Explainable machine learning predictions to help anesthesiologists prevent hypoxemia during surgery

Authors: Scott M. Lundberg¹, Bala Nair^{2,6}, Monica S. Vavilala^{2,6}, Mayumi Horibe⁴, Michael J. Eisses^{2,3}, Trevor Adams^{2,3}, David E. Liston^{2,3}, Daniel King-Wai Low^{2,3}, Shu-Fang Newman², Jerry Kim^{2,3}, Su-In Lee^{1,5}*

Affiliations:

¹Paul G. Allen School of Computer Science and Engineering, University of Washington, Seattle, WA, USA.

One Sentence Summary: We present a new machine learning based system called *Prescience* that provides interpretable real-time predictions to help anesthesiologists prevent hypoxemia during surgery.

Abstract: Hypoxemia causes serious patient harm, and while anesthesiologists strive to avoid hypoxemia during surgery, anesthesiologists are not reliably able to predict which patients will have intraoperative hypoxemia. Using minute by minute EMR data from fifty thousand surgeries we developed and tested a machine learning based system called *Prescience* that predicts real-time hypoxemia risk and presents *an explanation of factors* contributing to that risk during general anesthesia. Prescience improved anesthesiologists' performance when providing *interpretable hypoxemia risks with contributing factors*. The results suggest that if anesthesiologists currently anticipate 15% of events, then with Prescience assistance they could anticipate 30% of events or an estimated additional 2.4 million annually in the US, a large portion of which may be preventable because they are attributable to modifiable factors. The prediction explanations are broadly consistent with the literature and anesthesiologists' prior knowledge. Prescience can also improve clinical understanding of hypoxemia risk during anesthesia by providing general insights into the exact changes in risk induced by certain patient or procedure characteristics. Making predictions of complex medical machine learning models (such as Prescience) interpretable has broad applicability to other data-driven prediction tasks in medicine.

²Department of Anesthesiology and Pain Medicine, University of Washington, Seattle, WA, USA.

³Seattle Children's Hospital, Seattle, WA, USA.

⁴Veterans Affairs Puget Sound Health Care System, Seattle, WA, USA.

⁵Department of Genome Sciences University of Washington, Seattle, WA, USA.

⁶Harborview Injury Prevention and Research Center, Seattle, WA, USA

^{*}To whom correspondence should be addressed: suinlee@cs.washington.edu

Introduction

Over 200 million surgeries are performed worldwide every year, with 30 million in the United States alone (1). Though an integral part of healthcare, surgery and anesthesia pose considerable risk of complications and death. Studies have shown a perioperative mortality rate of 0.4 to 0.8% and a complication rate of 3 to 17%, just in industrialized countries (2, 3). Fortunately, half of these complications are preventable (2, 3). With increasing adoption of electronic medical record systems, high fidelity heterogeneous data are being captured during surgery and anesthesia care, yet the utilization of these data to improve patient safety and quality of care remains poor (4). There is untapped potential for data science to utilize perioperative data to positively impact surgical and anesthesia care (5). Recognizing this unmet need and leveraging recent advances in perioperative informatics, we present new data science methods to predict harmful physiological events and so inform anesthesiologists.

Hypoxemia or low arterial blood oxygen tension is an unwanted physiological condition and can cause serious patient harm during general anesthesia and surgery (6). Hypoxemia is associated with cardiac arrest, cardiac arrhythmias, postoperative infections and wound healing impairments, decreased cognitive function and delirium, and cerebral ischemia through a number of metabolic pathways (7). Despite the advent and use of pulse oximetry to continuously monitor blood oxygen saturation (SpO₂) during general and regional anesthesia, hypoxemia can neither be reliably predicted nor prevented (8). Real-time blood oxygen monitoring through pulse oximetry only allows anesthesiologists to take reactive actions to minimize the duration of hypoxemic episodes after occurrence. If hypoxemia can be predicted or anticipated before it occurs, then actions can be taken by anesthesiologists to proactively prevent hypoxemia and minimize patient harm.

Machine learning (ML) techniques use statistical methods to infer relationships between patient features and outcomes in large datasets, and have been successfully applied to predict adverse events in health care settings, such as sepsis, or patient deterioration in the intensive care unit (9–12). Yet ML techniques to predict adverse events such as hypoxemia in a considerably more critical and complex setting such as the operating room are currently lacking. Moreover, though previous complex ML approaches provide good prediction accuracy, their application in an actual clinical setting is limited because their predictions are difficult to interpret, and hence not actionable. Interpretable methods explain *why* a certain prediction was made for a patient, i.e., specific patient characteristics that led to the prediction. This lack of interpretability has thus far limited the use of powerful methods such as deep learning and ensemble models in medical decision support.

We present an ensemble model based machine learning method, *Prescience*, that predicts the near-term risk of hypoxemia during surgery and explains the patient and surgery specific factors that led to that risk (Figure 1). We believe this is an important step forward for machine learning in medicine because while machine learning models have significantly improved the ability to predict a patient's future condition (13), the inability to explain the predictions from accurate, complex models is a serious limitation. Understanding what drives a prediction is important for determining targeted interventions in a clinical setting. For this reason, machine learning methods employed in clinical applications avoid using complex, yet more accurate models and retreat to simpler interpretable (e.g., linear) models at the expense of poorer accuracy. We demonstrate how to retain interpretability, even when complex models such as nonparametric methods or deep learning are used, by developing a method to provide theoretically justified explanations of model predictions building on recent advances in model-agnostic prediction explanation methods (14–16). This allows these accurate, but traditionally hard to interpret, models to be used while still providing intuitive explanations of what led to a patient's predicted risk. Our ability to provide simple explanations of predictions from arbitrarily complex models eliminates the typical accuracy vs. interpretability tradeoff, thus allowing broader applicability of machine learning to medicine.

Prescience was trained to use standard operating room sensors to predict hypoxemic events in the near future and explain why an event is, or is not, likely to occur. It departs from the relatively few previous approaches to this problem in two important ways: First, unlike ElMoaqet et al. who used a linear autoregressive support vector machine on arterial oxygen saturation times series (9), and Tarassenko et al. who used Parzen windows to find outliers from five input patient measurement types (17), Prescience integrates a comprehensive dataset from a hospital's Anesthesia Information Management System (AIMS) (see Methods for details). The AIMS data consists of high fidelity *real-time data* – such as time series data from patient monitors and anesthesia machines, bolus and infusion medications, input and output fluid totals, laboratory results, templated and free text descriptions of anesthesia techniques and management, and *static data* – such as American Society for Anesthesiology (ASA) physical status, surgical procedure and diagnoses codes (18), as well as patient demographic information such as age, sex, smoking status, height and weight. By continuously integrating a broad set of patient and procedure *features* extracted from the AIMS data, Prescience surpasses human-level accuracy while maintaining consistent performance during every minute of a surgery.

Second, Prescience explains why a prediction was made, regardless of the complexity of the machine learning model used to make the prediction. Significant progress has been made recently integrating predictive machine learning solutions into medical care (9-12). However, accurately and intuitively conveying to doctors why a prediction was made remains a key challenge. For example, a numeric representation of risk (e.g., 2.4 odds ratio in Figure 1) is useful. However, a more detailed presentation that shows the risk is due to the patient's BMI, current tidal volume, and pulse rate is more clinically meaningful since some factors may be modifiable and result in clinical changes mitigating that risk.

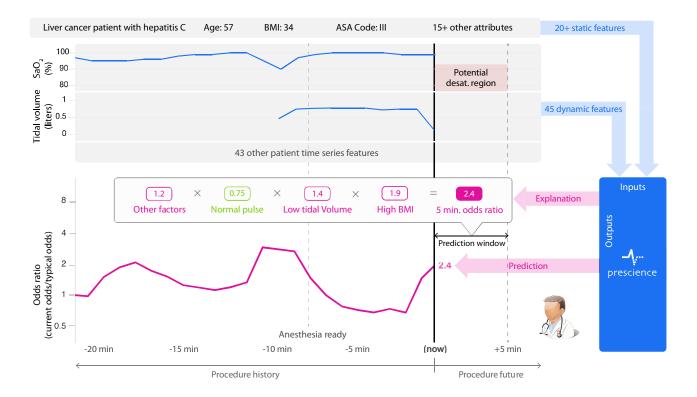


Fig 1. Prescience integrates many data sources into a single risk, which is explained through a succinct visual summary. A wide variety of data sources were used to build a predictive model of hypoxemia events. An explanation (shown above) is then built for each prediction. Purple features have values that increased risk, while green features decreased hypoxemia risk. The combination of impacts of all features is the predicted Prescience risk; in this case the odds are 2.4 times higher than normal.

Results

To demonstrate the value of Prescience's explained predictions and gain insight into factors affecting intraoperative hypoxemia, we present the following results: 1) a comparison of Prescience hypoxemia predictions against anesthesiologists' predictions with and without the aid of Prescience, 2) an example of how Prescience explains hypoxemia risk at a specific time-point during a surgical procedure, 3) a comparative summary of relevant AIMS data features for hypoxemia prediction chosen by Prescience and by anesthesiologists, and 4) a detailed presentation of key risk factors for hypoxemia identified by Prescience.

Prescience overview – data preparation, model learning and feature importance estimation

Based on World Health Organization recommendations and for the purposes of prediction, we defined hypoxemia as the drop in SpO₂, i.e., arterial blood oxygen saturation as measured by pulse oximetry, to 92% or lower (see Methods; Supplementary Figure 1). From the AIMS data, we extracted 3,797 *static features* for each patient and an expanded set of 3,905 *real-time* and *static features* for each time point during surgery (see Methods; Supplementary Table 1). We excluded cases (heart transplant, lung transplant, tracheostomy, and coronary artery bypass surgeries) in which SpO₂ and other hemodynamic parameters can be significantly affected by non-physiological measurements such as during cardiopulmonary bypass. All the experiments were performed after appropriate Institutional Review Board (IRB) approval (see Methods), with clinical data summarized in Table 1.

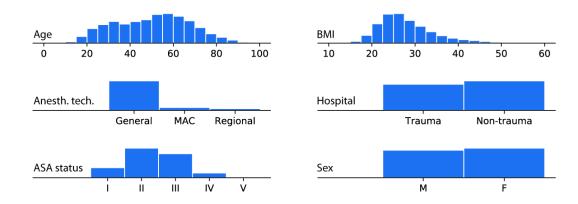


Table 1. Patient and procedure characteristics. Histograms summarizing basic properties of the anesthesia procedures used for training (ASA stands for American Society of Anesthesiologists, and BMI means Body Mass Index). Prescience was trained and evaluated using data from 53,126 procedures recorded at two hospitals over two years.

We trained a gradient boosting machine model (19) to solve the following two types of prediction problems: A) *initial prediction*: predicting at the start of a procedure the risk of hypoxemia anytime during a procedure based on the static features, and B) *real-time prediction*: predicting in the next 5 minutes at various points of the operative period based on real-time and static features collected up to that time point. For task A) we used 42,420 procedures as *training samples* to train the gradient boosting machine, 5,649 procedures as *validation samples* to choose the tuning parameters for the gradient boosting machine (and other prediction models for comparison), and 5,057 as *test samples* for comparing across different prediction models (Supplementary Figure 4). For task B), we used 8,087,476 time points as training samples, 1,053,629 as validation samples, and 963,674 as test samples, where all

time points from the same procedure were included in the same sample set (Supplementary Figure 3). To ensure that there was no bias towards the final test set, the test data was initially compressed and left compressed until method development was completed.

As shown in Supplementary Figures 3 and 4, the gradient boosting machine outperforms alternative prediction models previously used for similar problems.

For tasks A) and B), we use 198 and 523 test samples respectively for evaluating anesthesiologists' performance for initial and real-time prediction tasks, respectively (see below). Prescience outputs the risk prediction and its explanations (Figure 1; Figure 3A) which show a set of features that increased (purple-colored) and decreased (green-colored) the risk.

We developed an efficient, theoretically justified machine learning technique to estimate the importance of each feature on a prediction made for a single patient, which drives real-time explanations (Figure 3) for the Prescience model. We verified the quality of the explanations given to the anesthesiologists (in the experiments described below) by comparing the explanations with the change in model output when a feature is perturbed (Supplementary Figure 5). We also developed effective visualizations of these explanations that encodes them in a compact visual form for anesthesiologists (Figure 1; Supplementary Figures 6-8), and a more detailed visualization which highlights the relevant contributing features (Figure 3) (see Methods for details).

Prescience improves anesthesiologist's ability to predict hypoxemia

To test the potential of Prescience to aid hypoxemia prediction we replayed prerecorded surgery data from test sample procedures in a web-based visualization to five practicing anesthesiologists (Supplementary Figures 6-8). Each anesthesiologist was given two types of prediction tasks: A) *initial prediction* (198 tasks), and B) *real-time prediction* (523 tasks). For each prediction task, anesthesiologists were asked to provide a *relative risk* of hypoxemia as compared to a normal acceptable risk, for example, 0.01 for 1/100th the normal risk or 3.4 for 3.4 times the normal risk. These relative risks were then used to calculate standard receiver operating curves averaged over five anesthesiologists as shown in Figure 2 which plots the true positive rate (i.e., % of desaturations correctly predicted) in the y-axis against the false positive rate (i.e., % of non-desaturations incorrectly predicted) in the x-axis.

Figure 2A-B shows that for both types of prediction tasks, predictions made by Prescience (purple) are considerably more accurate than anesthesiologists' predictions (green). The prediction accuracy of anesthesiologists (green) markedly improved when the anesthesiologists were given Prescience's risk prediction and its explanations in addition to the original procedure data (blue) (Supplementary Figures 6-8). A clear separation between the performance of anesthesiologists with (blue) and without (green) the aid of Prescience is observed for both initial prediction (Figure 2A, P-value < 0.0001) and real-time prediction (Figure 2B, P-value < 0.0001). This suggests that Prescience can enhance anesthesiologists' assessment of future risk and their ability to proactively anticipate hypoxemia events. Interestingly, the prediction performance of anesthesiologists with Prescience explanations (blue) was slightly lower than direct predictions from Prescience (purple). This means that when the anesthesiologists adjust their risk estimate for a patient away from what Prescience originally predicted they are more likely to be wrong than right.

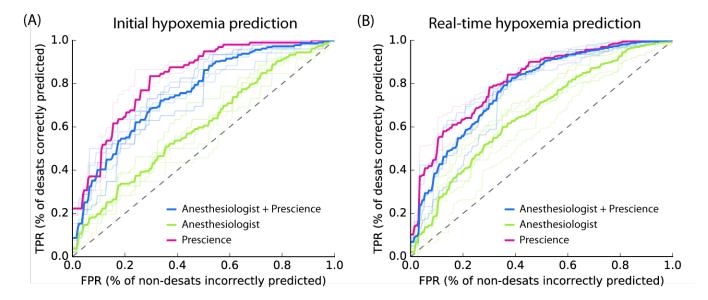


Fig 2. Pooled comparison of five anesthesiologists' prediction performance with and without assistance by Prescience. Receiver Operating Characteristic (ROC) plots comparing five anesthesiologists' predictions from recorded data with and without Prescience assistance. Light colored lines represent individual anesthesiologist's performances; dark lines represent their average performance. (**A**) For initial risk prediction, anesthesiologists (green, AUC = 0.60) performed significantly better with Prescience assistance (blue, AUC = 0.76; P-value < 0.0001) than without Prescience assistance, and Prescience performed better in a direct comparison with anesthesiologists (purple AUC = 0.83); P-value = 0.0003). (**B**) For intraoperative real-time (next 5 minute) risk prediction anesthesiologists (green, AUC = 0.66) again performed better with Prescience assistance (blue, AUC = 0.78; P-value < 0.0001), and Prescience alone outperformed anesthesiologists predictions (purple, AUC = 0.81); P-value = 0.0068). Note that The False Positive Rate (FPR) (x-axis) measures how many points without upcoming hypoxemia were incorrectly predicted to have upcoming hypoxemia. The True Positive Rate (TPR) (y-axis) measures what percentage of hypoxemic events was correctly predicted. P-values were computed using bootstrap resampling over the tested time points.

To avoid the scenario in which an anesthesiologist is tested on the same prediction task twice – one with and the other without Prescience, we created replicate test sets by dividing the prediction tasks into two groups of similar size: (100, 98) tasks for initial prediction and (260, 263) tasks for real-time prediction. Each of the five recruited anesthesiologists was assigned to receive Prescience's assistance in one of these two replicate test sets (Methods). The procedures shown to anesthesiologists were chosen such that ~50% showed at least one incident of hypoxemia (for preoperative prediction), and time points were chosen such that ~33% had hypoxemia in the next 5 minutes (for intraoperative prediction).

If we extrapolate the real-time results to the 30 million annual surgeries in the US under the assumption that doctors anticipate 15% of hypoxemic events while SpO_2 is still ≥ 95 , then with Prescience assistance they may be able to anticipate 30% of these events, or approximately 2.4 million additional episodes of hypoxemia annually (defined here as $SpO_2 \leq 92$). Since 20% percent of the Prescience risk prediction is based on drugs and settings under the control of the anesthesiologist (Supplementary Table 5; $Supp_RealtimeFeatureTable_5.csv$), a large portion of these predicted events may be preventable.

Explained risks reveal both procedure and time specific effects

An explanation from Prescience represents the effects of interpretable groups of patient features (see Figure 1 and Figure 3A). These effects explain why the model predicted a specific risk, and thus allow an anesthesiologist to plan appropriate interventions. In Figure 1 only the most significant features contributing to hypoxemia risk are shown for quick reference, however in Figure 3 the relative contributions of all patient and case features (i.e. attributes) towards hypoxemia risk can be seen at every sample time point during a procedure (Figure 3B). Without a meaningful explanation, the sudden increase in risk shown at the time point marked 'Now' might be hard to interpret; however, by representing the predicted risk as a cumulative effect of contributing patient and procedure features, the reason for the increase becomes clear (Figure 3A).

The increase in the risk of hypoxemia in the next 5 minutes shown in Figure 3 is driven by a set of features capturing both static attributes, such as patient height and weight, and dynamic parametric values, such as tidal volume (i.e., volume of gas exhaled per breath) and administration of drugs. The risk explanation bar in Figure 3A has purple-colored features that push the risk higher (to the right) and green-colored features that push the risk lower (to the left). Each group of features is sorted by the magnitude of their impact, and the largest impact features are labeled. Through this representation we can see that many of the 3,905 real-time features have only a small impact, and the risk for this time point is predominantly driven by a few features. Figure 3B shows the trend in the Prescience risk predictions over the course of the procedure. The plot in Figure 3B is equivalent to rotating the feature explanation in Figure 3A by 90 degrees and then stacking the explanations for each time point horizontally. We can see from the risk trend in Figure 3B that the large increase in risk at the current time was driven by 'Tidal volume', meaning a drop in the patient's tidal volume. The future SpO₂ (blood oxygen concentration) measurements confirm that the patient did indeed progress to hypoxemia (i.e., $SpO_2 \le 92$). Not only does Prescience alert anesthesiologists when a patient's risk for hypoxemia is high, but also provides information on the factors and their relative contributions driving the risk. This informed risk prediction enables anesthesiologists to plan an appropriate course of action to avoid hypoxemia.

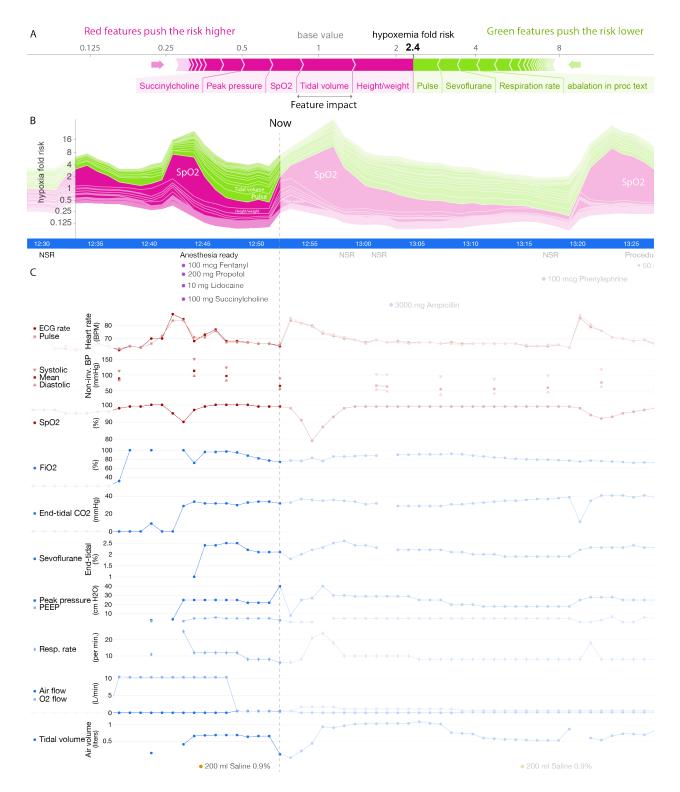


Fig 3. Sample real-time prediction during a procedure. One hour of data is shown from a procedure. (A) Explained risk of hypoxemia in the next five minutes. (B) Plot of the explained risks evolving over time. This plot is equivalent to rotating (A) 90 degrees and stacking the risk explanations for every time point horizontally. (C) A subset of the patient data for this procedure, plotted both before and after the current time point.

Averaged feature importance estimates broadly align with a survey of prior expectations

To gain an understanding of the general impact of features across all procedures we computed the average importance of each feature in the Prescience model. In contrast to the explanations shown in Figure 1 and Figure 3A which are specific to a single prediction at a particular time point, these average feature importance estimates are over many procedures and time points (20). These averaged feature importance estimates are shown for both initial prediction (Figure 4A) and real-time prediction (Figure 4B).

To estimate which clinical features anesthesiologists use to estimate hypoxemia risk, we first performed a survey before using Prescience, which asked four anesthesiologists to list the most important factors they consider when assessing the risk of hypoxemia, both before (for initial prediction) and during a procedure (for real-time prediction). Their responses were then aggregated into a single ranked list of features (Supplementary Tables 2 and 3). Figure 4 shows the rankings chosen by anesthesiologists next to the feature importance estimates derived by Prescience for (A) initial and (B) real-time predictions. The ranking of features by anesthesiologists appears to correspond well with the ranking by Prescience.

As another way to measure which features anesthesiologists' think contribute to hypoxemia, we learned from anesthesiologists' behavior by training a separate gradient boosting machine model based on their predictions. This allows a direct comparison between the anesthesiologists and Prescience on the same set of features. We fit this model to all the anesthesiologist relative risk predictions using 10-fold cross validation. We then computed the feature importance estimates for this model that was trained to mimic the behavior of anesthesiologists. Given the relatively smaller set of training examples used to train the model (198 initial predictions, and 523 real-time predictions), we used bootstrapping to estimate the variability of the feature importance estimates (Figure 4 right).


In general, there is reasonable agreement between the Prescience feature importance estimates and those identified by the anesthesiologists. However, there are important differences that may stem from the comprehensive nature of the Prescience analysis, while anesthesiologists necessarily focus on what they consider the most likely causes for hypoxemia concern. One striking difference is the reduced role of current SpO_2 levels in anesthesiologists' predictions. While anesthesiologists are clearly influenced by the recent patterns of patient SpO_2 levels, Prescience strongly depends on these patterns, while anesthesiologists appear to be equally influenced by other factors, such as end tidal CO_2 and peak ventilation pressure.

Our study used data from two hospitals and initial hypoxemia predictions were driven by a bias between the two hospitals. This is perhaps unsurprising since one hospital is a level 1 trauma center and a significant proportion of its surgical cases involve trauma patients who are more susceptible to hypoxemia. However, it is interesting to note that the importance of hospital as a risk factor became insignificant for the intraoperative real-time predictions, presumably because the risk differences in each hospital were captured by the real-time features.

Among the static features, BMI (body mass index) and age were significant risk factors. These features are well understood in the medical literature as risk factors that can increase the chances of hypoxemia (21, 22). The American Society of Anesthesiology (ASA) physical status feature represents the severity of a patient's medical condition and a higher ASA number indicates a higher comorbidity. Prescience determined that higher ASA physical status values predisposes a patient to higher hypoxemia risk. While this finding may be clinically intuitive, anesthesiologists can now use this information in their preoperative evaluation as a pre-specified risk factor for intraoperative hypoxemia. Eye procedures were informative to the model and carried a reduced risk of hypoxemia, while surgeries for fractures had a slightly higher risk. These patterns may reflect the composite risk of hypoxemia to patients undergoing these particular procedures. In the case of eye surgeries, the risk was lower even though many are elderly and have accompanying co-morbid conditions. The low risk of eye procedures is in contrast to

fractures which carry an increased risk for hypoxemia. Together these findings provide new data on the relative "risk" of these procedures which has implications for anesthesia staffing, need for equipment, and preparation for the ability to rescue patients from hypoxemia. Eye procedures and surgeries for fractures are two examples of text based features extracted from diagnosis and preoperative procedure notes. They demonstrate that unstructured text notes can be combined with structured patient data to improve patient risk prediction. Although many of the risk factors identified by Prescience reconfirmed those expected by the anesthesiologists, it is informative that Prescience independently identified these features with no prior knowledge.

Among real-time (intraoperative) features SpO₂ (arterial oxygen saturation) is, as expected, the strongest predictor of future potential drops in SpO₂. End tidal CO₂ (amount of carbon dioxide exhaled by the patient) was also a significant intraoperative feature identified by Prescience as predictive of hypoxemia. Lower values may indicate inadequate ventilation or airway obstruction which can turn increase the risk of hypoxemia. Prescience also determined that hypotension (systolic blood pressure below 80) increases the risk of hypoxemia. On the other hand, increased FiO₂ (inspired O₂ concentration) and adequate positive pressure ventilation can reduce the risk of hypoxemia, as expected by the anesthesiologists.

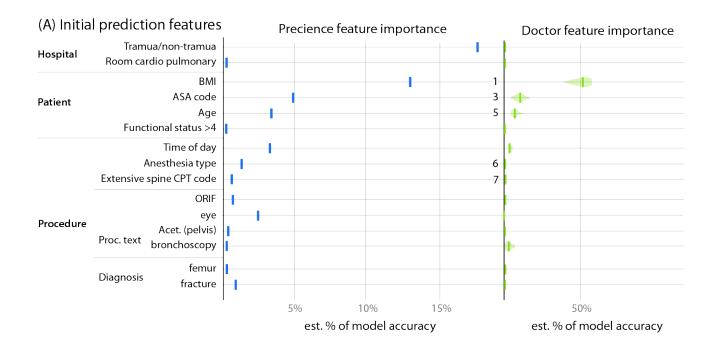

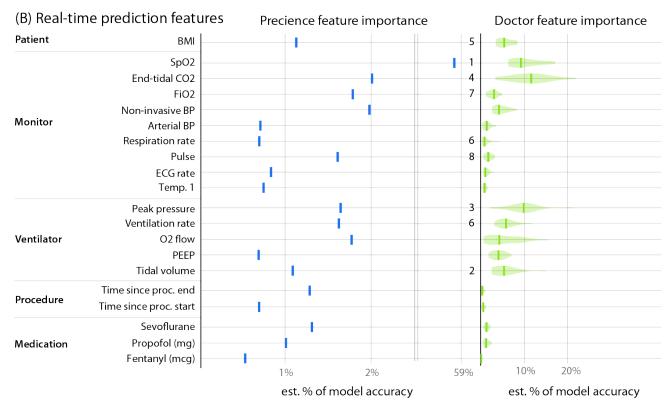


Fig 4. Comparison of averaged feature importance estimates between Prescience and anesthesiologists for both initial and real-time prediction. Importance estimates assigned by the Prescience model and anesthesiologists to the top features in both (A) initial and (B) real-time prediction. The importance of features is measured as the estimated percent of the model's prediction accuracy that is due to that feature. The numbers presented to the left of the imputed anesthesiologist importance estimates are feature rankings from a consensus of anesthesiologist responses about which features they believed would be important. The quantitative

anesthesiologist feature importance estimates were estimated using 20 bootstrapped models trained to mimic the anesthesiologist's predictions when unassisted by Prescience. ¹

Prescience's estimated importance of individual features on hypoxemia risk highlight important clinical relationships

Three important features each for beginning of surgery (initial) and during surgery (real-time) predictions were chosen to illustrate how the Prescience model modifies hypoxemia risk based on changes to feature characteristics (Figure 4). While many such relationships are present for the various features, Figure 5 shows a representative selection demonstrating informative risk relationships that are captured in the Prescience model.

Among static features we find that patient BMI has a clear effect on the risk of hypoxemia. When the BMI is over 26, the risk of hypoxemia increases linearly until it has more than doubled when BMI is over 50. Though a qualitative association between hypoxemia and body weight is well established in the field of anesthesia (21, 22), Prescience quantifies this relative risk.

Prescience shows that patients with higher ASA physical status codes have higher risk of intraoperative hypoxemia. This is not surprising since higher ASA codes represent increased severity of a patent's physical condition such as preexisting pulmonary and cardiac conditions that can predispose a patient to develop hypoxemia. Prescience data support clinical observations that the risk of hypoxemia more than doubles when the ASA code increases from I to V. Advancing age also predicted intraoperative hypoxemia, likely representing the presence of comorbidities (21). These data show that BMI > 30, which meets the clinical definition of obesity (23), is associated with intraoperative hypoxemia, suggesting impaired pulmonary mechanics. While we agree that these findings confirm clinical observations and suspicions of the relationship between these patient factors and adverse anesthesiology outcomes, Prescience quantifies this association and the risks, giving a more clinically useful interpretation to anesthesiologists.

For real-time prediction, measurements from each time series are represented by a set of multiple features. For simplicity, we focus here only on the effect of the shortest time lag exponentially weighted moving average, which essentially represents the most recent reported value in the time series (see Methods for details).

Tidal volume represents the amount of gas exhaled per breath when the patient is either breathing spontaneously or mechanically ventilated during general anesthesia. As the tidal volume drops below 0.6 liters (keeping all other features the same), Prescience risk for hypoxemia increases. This increase could be due to hypoventilation, in which case anesthesiologists take preventative steps to avoid inadequate ventilation.

End tidal CO₂ represents the amount of carbon dioxide exhaled gas. Figure 5 shows the relationship between end tidal CO₂ and risk of hypoxemia under general anesthesia. End tidal CO₂ below 35 mmHg is associated with an increasing risk of intraoperative hypoxemia. While we cannot definitively attribute hypoxapnia with intraoperative hypoxemia, these associations may represent underlying patient conditions such as chronic obstructive pulmonary disease that affect both physiological conditions. Alternately, the low end tidal CO₂ and may result from either intentional or unwanted hyperventilation during anesthesia care.

¹ In Figure 4A the 2nd and 4th features chosen by anesthesiologists (lung disease and asthma, respectively) were not well captured in the AIMS data, and did not show up as important in the Prescience model.

Examining FiO_2 is important because anesthesiologists can control the amount of oxygen delivered to patients. Current practice is to not provide all patients with 100% FiO_2 since not all patients need it, prolonged ventilation with 100% FiO_2 is associated with pulmonary atelectasis, and because oxygen if delivered when not needed is costly and wasteful. These data show that FiO_2 below 40% is independently associated with intraoperative hypoxemia irrespective of other features. These findings provide important information regarding safe practice of FiO_2 in patients during general anesthesia. It is possible that the routine practice of maintaining FiO_2 30% or close to room air may be harmful to patients and not desirable. While these effects are adjusted for all other available features, it is important to note that as with any observational study some residual confounding with patient risk may still exist. This could explain the increase in hypoxemia risk we observed for high O_2 levels.

These representative features illustrate the power of our machine learning-based prediction method, Prescience, to not only provide explained risk predictions, but also quantitative insights into the exact change in risk induced by certain patient or procedure characteristics.

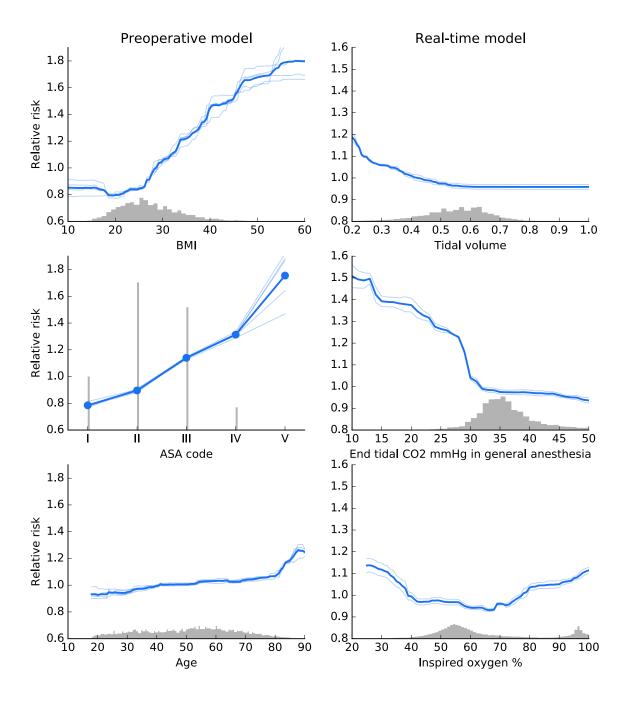


Fig 5. Effect of varying individual feature values for both initial features (left) and real-time features (right). These partial dependence plots show the change in hypoxemia risk for all values of a given feature. The gray histograms on each plot show the distribution of values for that feature in the validation dataset. Light colored lines represent model variability from bootstrap resampling of the training data.

Discussion

To the best of our knowledge, Prescience is the first method designed to comprehensively integrate high fidelity operating room data to predict intraoperative hypoxemia events before they occur. Based on a comparison against practicing anesthesiologists and existing computational methods applied to other clinical problems, Prescience achieves superior performance when predicting hypoxemia risk from electronically recorded data.

To our knowledge, Prescience is also the first method to combine high accuracy complex models with interpretable explanations. This combination of accuracy and interpretability allows physicians to receive the best possible predictions while also gaining insight into why those predictions were made. To test how Prescience predictions with explanations would impact an anesthesiologist's ability to estimate hypoxemia risk we compared anesthesiologist predictions with and without Prescience assistance. We observed a clear increase in prediction accuracy when doctors were assisted by Prescience, demonstrating that anesthesiologists may make more accurate hypoxemia risk assessments in the operating room if they had access to Prescience.

It should be clarified that our exercise at developing machine learning methods to predict intraoperative hypoxemia, though promising, should still be considered an initial attempt. In this first attempt, we did not categorize procedures to assess hypoxemia predictions in specific types of procedures. For this reason, clinical interpretation of the results had to be somewhat generic. For enhanced interpretation of risks, future attempts need to focus on specific categories of cases and phases of anesthesia. Another future enhancement would be integrating additional preoperative data such as patient's detailed medical history into the prediction models. Higher fidelity intraoperative data such as patient monitor waveform data could enrich machine learning thus potentially leading to more accurate predictions. Prospective trials of Prescience during live procedures are also needed to verify the improvements in anesthesiologist's performance we retrospectively observed in prerecorded procedures.

This paper focuses on hypoxemia risk during intraoperative anesthesia care. However, the importance of coupling accurate predictions from complex models with interpretable explanations of why a prediction was made, has broad applicability throughout medicine. Extending the approach taken by Prescience and providing these model explanation tools to the community is a clear next step. Because Prescience effectively decouples the interpretable explanation from the prediction model, we are also free to continue to refine the core prediction model without changing the user experience for anesthesiologists.

The global risk profiles learned by Prescience (Figures 4-5) are clinically relevant for a number of reasons. First, they show that in the health system examined, trauma hospital patients may be more critically ill as they have more intraoperative hypoxemia. In current times when harmonization of care and standardization are considered to reduce unwanted clinical variation, these data suggest that resources may need to be differentially deployed to address differential rates of adverse events. Second, anesthesiologists can now quantify risks of intraoperative hypoxemia adjusted for other factors to the very elderly, those who are overweight, and those with more comorbid conditions. The exact relationships described in Figure 5 clearly show the patterns and threshold points for the risk. Whereas low tidal volume is often suggested for patients with acute lung injury (24), these data suggest that overall, low lung tidal volumes are, in-fact, associated with intraoperative hypoxemia. The relationship between low end-tidal CO₂ levels and intraoperative hypoxemia may reflect underlying critical illness. Despite our inability to fully exclude residual confounding, these data shed new light on physiological relationships as well as provide a mechanism to facilitate provision of anesthesia care that can mitigate intraoperative hypoxemia.

The field of medicine is full of many exciting data science challenges that have the potential to fundamentally impact the way medicine is practiced. More and more data driven predictions of patient outcomes are being proposed and used. However, black-box prediction models which provide simply predictions, without explanation, are difficult for physicians to trust and provide little insight about how they should respond. The interpretable explanations used by Prescience represent a powerful technique that can transform any current prediction method from one that provides *what* the prediction is, into one that also explains *why*.

Materials and Methods

IRB statement

The electronic data for this study was retrieved from institutional electronic medical record and data warehouse systems after receiving approval from the Institutional Review Board (University of Washington Human Subjects Division, Approval #46889). Protected health information was excluded from the data set that was used for machine learning methods.

Data sources

Our hospital system has installed an Anesthesia Information Management System (AIMS) (Merge AIM, Merge Inc, Hartland, WI) that automatically captures minute by minute hemodynamic and ventilation parameters from the patient monitor and the anesthesia machine. The system also integrates with other hospital electronic medical record (EMR) systems to automatically acquire laboratory and patient registration information. The automatic capture of data is supplemented by manual documentation of medications and anesthesia interventions to complete the anesthesia record during a surgical episode. For the current project, we extracted the high-fidelity anesthesia data from the AIMS database for the period May 2012 through June 2014. Additionally, for each patient, medical history data were extracted from our EMR data warehouse (Caradigm, Bellevue, WA). The high-fidelity anesthesia record data and the corresponding medical history data from the hospital EMR formed the underlying data for machine learning. The various data elements used for machine learning are outlined in Supplementary Table 1.

Patient Age Age of patient Sex Sex ASA Physical Status ASA physical status Height (cm) Height of patient Weight (kg) Weight of patient Procedure Procedure Procedure Procedure of Procedure description Procedure code Billing codes Anesthesia Crosswalk/ procedure codes Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis odes (CD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Etant Time of Frocedure End Emergence Time of Procedure End Emergence Leave OR Anesthesia Td Time of leave OR/ Transport to recovery Anesthesia End Patient monitor/Ventilator Data NIBP—Sys Cuff BP—Systolic (Patient Monitor) NIBP—Sys Cuff BP—Systolic (Patient monitor) NIBP—Sys Cuff BF—Systolic (Patient monitor)	Element	Description
Sex Sex ASA Physical Status ASA physical status Height (cm) Height of patient Weight (kg) Weight of patient Patient class Inpatient or outpatient Procedure Procedure Procedure code Surgical procedure Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes ICD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start In Room Time of Anesthesia Start In Room Time of patient in room Induction Time of induction start Anesthesia Ready Time of induction start Anesthesia Ready Time of procedure start (incision) Closing Time of Procedure End Emergence Time of leave OR/ Transport to recovery A	Patient	
ASA physical status Height (cm) Height of patient Weight (kg) Weight of patient Procedure Procedure Procedure Procedure Procedure Surgical procedure Billing codes Diagnosis Diagnosis description Posality Hospital facility Emergency status Anesthesia type Case Events In Room Induction In Room Induction In Time of Induction end or Anesthesia ready Procedure Start Closing Time of Procedure start (incision) Closing Procedure Start Im of Procedure tart of emergence Leave OR Ime of Anesthesia End Patient Monitor/Ventilator Data Heart Rate Quiser in Monitor (Patient Monitor) Pulse rate from pulse oximetry (Patient Monitor) Pulse rate	Age	Age of patient
Height (cm) Height of patient Weight (kg) Weight of patient Protedure Procedure Procedure Procedure Procedure Surgical procedure Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes ICD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start In mo of Anesthesia Start In Room Time of Induction start Anesthesia Ready Time of Procedure and vanishesia ready Procedure Start Time of Procedure start (incision) Closing Time of Procedure End Emergence In Time of Sart of emergence Leave OR Time of leave OR, Transport to recovery Anesthesia End Patient monitor/Ventilator Data Pulse rate from EUG signal (Patient monitor) Og. Sat Pulse rate from pulse oximetry (Patient Monitor) Pulse rate from pulse oximetry (Patient Monitor)		<u> </u>
Weight (kg) Weight of patient Patient class Inpatient or outpatient Procedure Procedure code Procedure description Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Procedure start (incision) Closing Time of Frocedure End Emergence Time of start of emergence Leave OR Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of mea	ASA Physical Status	ASA physical status
Patient class Inpatient or outpatient Procedure Procedure Procedure Procedure Procedure description Procedure code Surgical procedure Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes ICD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of start of emergence Leave OR Time of Leave OR/ Transport to recovery Anesthesia End Patient monitor/Ventilator Data Pulse rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Procedure Monitor) Pulse rate from pulse oximetry (Patient Monitor)	Height (cm)	Height of patient
Procedure Procedure Procedure description Procedure code Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis des ICD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Procedure End Emergence Time of start of emergence Leave OR Time of start of emergence Leave OR Time of Anesthesia End Patient monitor/Ventilator Data Heart Rate Heart rate from ECG signal (Patient monitor) Pulse rate from pulse oximetry (Patient Monitor) Pulse rate from pulse oximetry (Patient Monitor)	Weight (kg)	Weight of patient
Procedure Procedure description Procedure code Surgical procedure Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Ime of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Procedure start (incision) Closing Time of Procedure start (incision) Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of Anesthesia End Patient monitor/Ventilator Data Pulse rate from ECG signal (Patient Monitor) Pulse rate Pulse rate Pulse rate from pulse oximetry (Patient Monitor) Pulse rate	Patient class	Inpatient or outpatient
Procedure code Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes Location Operating room location Facility Hospital facility Emergency status Anesthesia type Type of anesthesia Anesthesia Start In Room Induction Induction Induction Induction Induction Induction Induction Induction In Forecedure Start In Forecedure Start In Forecedure Start In Fore of Induction end or Anesthesia ready Procedure Start Inime of Procedure start (incision) Closing Time of Closing Time of Closing Time of Procedure End Emergence Leave OR Anesthesia End Patient monitor/Ventilator Data Value, Time & Unit of measurement) Feat Rate Pulse rate Pulse rate from pulse oximetry (Patient Monitor) Pulse rate	Procedure	
Billing codes Anesthesia Crosswalk/ procedure codes Diagnosis Diagnosis description Diagnosis codes ICD-9/10 codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Closing Procedure Start Time of Procedure start (incision) Closing Time of Procedure End Emergence Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Procedure	Procedure description
Diagnosis Diagnosis description Diagnosis codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Ime of Anesthesia Start Ime of patient in room Induction Ime of Induction start Anesthesia Ready Ime of Induction end or Anesthesia ready Procedure Start Ime of Procedure start (incision) Closing Ime of Start of emergence Emergence Ime of start of emergence Leave OR Ime of leave OR/ Transport to recovery Anesthesia End Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor) Pulse rate	Procedure code	Surgical procedure
Diagnosis codes Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start In Room Induction Induct	Billing codes	Anesthesia Crosswalk/ procedure codes
Location Operating room location Facility Hospital facility Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Time of Procedure End Emergence Time of Start of emergence Leave OR Time of Anesthesia End Patient monitor/Ventilator Data Heart Rate Heart rate from ECG signal (Patient Monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Diagnosis	Diagnosis description
Facility Emergency status Anesthesia type Type of anesthesia Case Events Anesthesia Start In Room Induction Induction Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Procedure End Emergence Leave OR Anesthesia End Patient monitor/Ventilator Data Heart Rate Pulse rate Vhether are from ECG signal (Patient Monitor) Pulse rate Pulse rate Vhether is an emergency or not – Y/N Anesthesia a nemergency or not – Y/N Type of anesthesia Time of Anesthesia Start Time of Anesthesia Start Time of Induction end or Anesthesia ready Time of Procedure start (incision) Time of Procedure End Time of Procedure End Time of Start of emergence Leave OR Time of start of emergence Leave OR Time of Anesthesia End Patient monitor/Ventilator Data Value, Time & Unit of measurement) Pulse rate from pulse oximetry (Patient Monitor) Pulse rate from pulse oximetry (Patient Monitor)	Diagnosis codes	ICD-9/10 codes
Emergency status Whether case is an emergency or not – Y/N Anesthesia type Type of anesthesia Case Events Anesthesia Start Time of Anesthesia Start In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Location	Operating room location
Anesthesia type Case Events Anesthesia Start In Room Induction Induction Induction In Forcedure Start In me of Induction end or Anesthesia ready Procedure Start In me of Procedure start (incision) Closing Induction Induction Induction Induction end or Anesthesia ready Procedure Start In me of Procedure start (incision) Induction Induction end or Anesthesia ready Procedure Start In me of Procedure start (incision) Induction end or Anesthesia ready Procedure Start In me of Procedure start (incision) Induction end or Anesthesia ready Procedure End In me of Closing In end Frocedure End In end	Facility	Hospital facility
Case Events Anesthesia Start In Room Induction end or Anesthesia ready Induction end or Anesthesia	Emergency status	Whether case is an emergency or not – Y/N
Anesthesia Start In Room Induction Time of patient in room Induction Anesthesia Ready Procedure Start Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Leave OR Time of start of emergence Leave OR Time of Anesthesia End Patient monitor/Ventilator Data Heart Rate Heart rate from ECG signal (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor) Pulse rate	Anesthesia type	Type of anesthesia
In Room Time of patient in room Induction Time of Induction start Anesthesia Ready Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Case Events	
Induction Anesthesia Ready Procedure Start Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Leave OR Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Anesthesia Start	Time of Anesthesia Start
Anesthesia Ready Procedure Start Time of Induction end or Anesthesia ready Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	In Room	Time of patient in room
Procedure Start Time of Procedure start (incision) Closing Time of Closing Procedure End Time of Procedure End Time of Procedure End Emergence Leave OR Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Induction	Time of Induction start
Closing Time of Closing Procedure End Time of Procedure End Emergence Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Anesthesia Ready	Time of Induction end or Anesthesia ready
Procedure End Emergence Leave OR Anesthesia End Time of Frocedure End Time of start of emergence Leave OR Time of leave OR/ Transport to recovery Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate	Procedure Start	Time of Procedure start (incision)
Emergence Leave OR Time of start of emergence Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Closing	Time of Closing
Leave OR Time of leave OR/ Transport to recovery Anesthesia End Time of Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate from pulse oximetry (Patient Monitor)	Procedure End	Time of Procedure End
Anesthesia End Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Emergence	Time of start of emergence
Patient monitor/Ventilator Data (Value, Time & Unit of measurement) Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Leave OR	. ,
Heart Rate Heart rate from ECG signal (Patient monitor) O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Anesthesia End	Time of Anesthesia End
O ₂ Sat O ₂ saturation from pulse oximetry (Patient Monitor) Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Patient monitor/Ventilator Data	(Value, Time & Unit of measurement)
Pulse rate Pulse rate from pulse oximetry (Patient Monitor)	Heart Rate	Heart rate from ECG signal (Patient monitor)
• • • • • • • • • • • • • • • • • • • •	O ₂ Sat	O ₂ saturation from pulse oximetry (Patient Monitor)
NIBP – Sys Cuff BP – systolic (Patient monitor)	Pulse rate	Pulse rate from pulse oximetry (Patient Monitor)
	NIBP – Sys	Cuff BP – systolic (Patient monitor)


NIBP - Dia Art BP - Sys Arterial BP - systolic (Patient monitor) Art BP - Dia Art BP - Dia Art BP - Mean Arterial BP - dasylic (Patient monitor) Art BP - Mean Arterial BP - dasylic (Patient monitor) Art BP - Mean Arterial BP - dasylic (Patient monitor) Art BP - Mean Arterial BP - dasylic (Patient monitor) Art BP - Mean Arterial BP - mean (Patient monitor) Art BP - Mean Arterial BP - mean (Patient monitor) PA - Sys Pulmonary artery pressure - dasylic (Patient monitor) CVP Central venous line pressure - mean (Patient monitor) EtCO, End tidal CO, (Capnography) (Patient monitor) EtCO, End tidal CO, (Capnography) (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Desflurane anesthetic agent (Patient monitor) ET N20 End tidal Soflurane anesthetic agent (Patient monitor) ET N20 End tidal Soflurane anesthetic agent (Patient monitor) BIS Bispectral Index (Patient Miss monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) FMP Temperature (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW Air flow rate (Ventilator) N,O FLOW O, flow rate (Ventilator) Type of cardiac rhythm Medications Delivery information of medications Delivery time (start / end for infusion medications) Drug Asme Dose Drug Asme Dose Drug Asme Drug Asse Fluid totals Time Time of fluid input or output Fluid totals Time Time of taking sample or Lab result time Lab result Lab result Lab result Lab result Unit Unit of measurement Notes Attestations, AlMS Clinical notes, Note option selections Time Time associated with the note	NURD Dia	Coff DD disabelia (Dationt associated)		
Art BP - Sys Arterial BP - systolic (Patient monitor) Art BP - Dia Art BP - Dia Arterial BP - diastolic (Patient monitor) Art BP - Mean Arterial BP - diastolic (Patient monitor) PA - Sys Pulmonary artery pressure - systolic (Patient monitor) PA - Dia PA - Dia Pulmonary artery pressure - diastolic (Patient monitor) CVP Central venous line pressure - mean (Patient monitor) EtCO; End tidal CO; (Capnography) (Patient monitor) EtCO; End tidal CO; (Capnography) (Patient monitor) ET CO; Inspired O; (Patient monitor) ET Sevo End tidal Sevoffurance anesthetic agent (Patient monitor) ET Des End tidal Sevoffurance anesthetic agent (Patient monitor) ET N20 End tidal Informance anesthetic agent (Patient monitor) ET N20 End tidal Nitrous oxide (Patient monitor) BIS Bis Bispectral Index (Patient/BIS monitor) SOI Signal quality index of BIS (Patient/BIS monitor) SOQ; Mixed venous oxygenation (Patient monitor) TEMP Temperature (Patient monitor) TOY Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PEP Peak Inspiratory pressure (Ventilator) PEP Peak Inspiratory pressure (Ventilator) PEP Peak Inspiratory pressure (Ventilator) No, O FLOW Air FLOW Air FLOW Air FLOW Air Flow at (Ventilator) No, O FLOW Air Flow are (Ventilator) No, O Flow are (Ventilator) No, O Flow are (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Positive End Expiration Pressure (Ventilator) Preg Positive End Expiration Pressure (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Peak Inspiratory pressure (Ventilator) Preg Positive End Expiration Pressure (Ventilator) Preg Positive End E		,		
Art BP - Dia Art BP - Mean Arterial BP - mean (Patient monitor) Art BP - Mean Arterial BP - mean (Patient monitor) PA - Sys Pulmonary artery pressure - systolic (Patient monitor) PA - Dia Pulmonary artery pressure - systolic (Patient monitor) PA - Dia Pulmonary artery pressure - diastolic (Patient monitor) ECCO End tidal CO; (Capnography) (Patient monitor) ECCO, End tidal CO; (Capnography) (Patient monitor) ET Sevo End tidal Sevidurane anesthetic agent (Patient monitor) ET Sevo End tidal Sevidurane anesthetic agent (Patient monitor) ET Sevo End tidal Desflurane anesthetic agent (Patient monitor) ET ISO End tidal Desflurane anesthetic agent (Patient monitor) ET ISO End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Introus oxide (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) ET N20 End tidal Sevidurane anesthetic agent (Patient monitor) EV Temperature (Patient monitor) EV Temperature (Patient monitor) EV Temperature (Patient monitor) EV Temperature (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) PEEP Postitue End Expiration Pressure (Ventilator) PEEP Postitue End Expiration Pressure (Ventilator) Op ELOW Air FLOW Air FLOW Air Flow rate (Ventilator) Air FLOW Air Flow rate (Ventilator) Air Flow rate (Ventilator) PEEP Postitue End Expiration Pressure (Ventilator) PEEP Postitue End Expiration Pressure (Ventilator) Porug Name Pose Drug Name Porug Nam		, ,		
Art BP - Mean	•	, , , , , , , , , , , , , , , , , , , ,		
PA - Sys Pulmonary artery pressure - systolic (Patient monitor) PA - Dia Pulmonary artery pressure - diastolic (Patient monitor) CVP Central venous line pressure - mean (Patient monitor) EtCO; End tidal CO; (Capnography) (Patient monitor) Resp Rate Measured respiration rate - capnography (Patient monitor) FIO; Inspired O; (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Sevoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET N2O End tidal Isoflurane anesthetic agent (Patient monitor) BIS BiS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BiS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO; Mixed venous oxygenation (Patient monitor) Vor Tidal volume (Patient monitor) Vor Tidal volume (Patient monitor) Vor Tidal volume (Patient monitor) Vor RATE Ventilator rate setting (Ventilator) PEP Positive End Expiration Pressure (Ventilator) PEP Positive End Expiration Pressure (Ventilator) Air FLOW O; flow rate (Ventilator) Air FLOW O; flow rate (Ventilator) N; O FLOW O; flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Time Delivery time (start / end for infusion medications) Drug Name Dose Dose Dose Drug dose Dose Unit Drug dose unit Put dotals Time Time of fluid input or output Fluid Name Fluid Name Fluid volume Fluid volume Fluid volume recorded Laboratory results Intraoperative laboratory results Intraoperative laboratory results Intraoperative laboratory results Intraoperative laboratory results Lab result Lab result Unit Unit of measurement Notes		, ,		
PA - Dia Pulmonary artery pressure - diastolic (Patient monitor) CVP Central venous line pressure - mean (Patient monitor) EtCO; End tidal CO; (Capnography) (Patient monitor) Resp Bate Measured respiration rate - capnography (Patient monitor) FIO; Inspired O; (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Desflurane anesthetic agent (Patient monitor) ET NZO End tidal Isoflurane anesthetic agent (Patient monitor) ET NZO End tidal Isoflurane anesthetic agent (Patient monitor) BIS Bispectral index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO ₂ Misked venous oxygenation (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) PIP Peak inspiratory pressure (Ventilator) PIP Peak inspiratory pressure (Ventilator) O; FLOW O; flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery time (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Drug Name Dose Drug Name Dose Drug dose Drug dose Drug dose Drug dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Time Time of fluid input or output Fluid Name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Ime of taking sample or Lab result time Lab result Unit Unit of measurement Notes Attestations, AlMS Clinical notes, Note option selections		,		
CVP EtCO; End tidal CO; (Capnography) (Patient monitor) EtCO; End tidal CO; (Capnography) (Patient monitor) Resp Rate Measured respiration rate - capnography (Patient monitor) FiO; Inspired O; (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET N2O End tidal Isoflurane anesthetic agent (Patient monitor) BIS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) FEMP Temperature (Patient monitor) FEMP Temperature (Patient monitor) TV Tidal volume (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIPP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW O, flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N,O FLOW O, flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N,O FLOW N,O Flow rate (Ventilator) Delivery information of medications Time Delivery ime (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Time Drug Asme Dose Drug dose Drug dose Drug dose Drug dose Drug dose Drug dose Drug dose Drug dose Drug dose Drug dose Lab result totals Time Time of fluid input or output Fluid Name Fluid Name Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result Result value Unit Unit of measurement Notes				
EtCO2 Resp Rate Measured respiration rate — capnography (Patient monitor) Resp Rate Measured respiration rate — capnography (Patient monitor) FIO2 Inspired O2, (Patient monitor) ET Sevo End tidal Sevoffurane anesthetic agent (Patient monitor) ET Des End tidal Sevoffurane anesthetic agent (Patient monitor) ET Des End tidal Instruane anesthetic agent (Patient monitor) ET N20 End tidal Instruane anesthetic agent (Patient monitor) BIS Bispectral Index (Patient Monitor) BIS Bispectral Index (Patient Monitor) SQI Signal quality index of BIS (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) FEMP Temperature (Patient monitor) FEMP Temperature (Patient monitor) CCO Continuous cardiac output (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N.O FLOW N.O Flow rate (Ventilator) POR O2 flow rate (Ventilator) Delivery information of medications Time Delivery time (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Time Delivery time (Start / end for infusion medications) Time Delivery Name Dose Drug Aose Dose Unit Drug Aose Dose Unit Drug Aose Drug Aose unit Drug dose unit Drug Aose		Pulmonary artery pressure – diastolic (Patient monitor)		
Resp Rate FiG.2 Inspired O. (Patient monitor) FiG.2 Inspired O. (Patient monitor) ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Desflurane anesthetic agent (Patient monitor) ET ISO End tidal Inspired O. (Patient monitor) ET ISO End tidal Inspired O. End tidal Inspired O. (Patient monitor) ET ISO End tidal Inspired anesthetic agent (Patient monitor) ET IXO End tidal Inspired anesthetic agent (Patient monitor) BIS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) CCC Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) PIP Peak Inspiratory pressure (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Opositive End Expiration Pressure (Ventilator) O.; FLOW O, Goverate (Ventilator) Air FLOW Air flow rate (Ventilator) N.; O FLOW Air (Ventilator) N.; O FLOW N.; O Flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Pedications Drug Name Delivery information of medications Time Delivery information of medications Time Delivery information of medications Time Dose Drug dose Drug dose Drug dose Dose Unit Drug dose unit Pour Gose Drug dose Dose Drug dose Dose Drug dose unit Fluid totals Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid volume recorded Laboratory results Lab result Lab result time Lab result Lab result time Lab result Result value Unit Unit Onessurement Notes		· · · · · · · · · · · · · · · · · · ·		
FIO2 ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Sevoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) BIS BIS Bispectral Index (Patient!/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PIP Pes Positive End Expiration Pressure (Ventilator) Air FLOW Air flow ate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Tye of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug Name Drug dose Drug dose Drug dose Drug dose Drug dose Unit of measurement Notes Notes Attestations, AIMS Clinical notes, Note option selections	EtCO ₂	End tidal CO₂ (Capnography) (Patient monitor)		
ET Sevo End tidal Sevoflurane anesthetic agent (Patient monitor) ET Des End tidal Desflurane anesthetic agent (Patient monitor) ET ISO End tidal Solfurane anesthetic agent (Patient monitor) ET N2O End tidal Nitrous oxide (Patient monitor) BIS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SvO2 Mixed venous oxygenation (Patient monitor) TCCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) TV V Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PPP Peak Inspiratory pressure (Ventilator) PPP Peak Inspiratory pressure (Ventilator) Air FLOW O, flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N2O FLOW O, Slow rate (Ventilator) N2O FLOW O, Flow rate (Ventilator) N2O FLOW N2O Flow rate (Ventilator) Time Delivery time (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Time Drug dose Dose Drug dose Dose Unit Drug dose Dose Unit Drug dose Dose Unit Drug dose Dose Unit Drug dose Drug dose Dose Unit Drug dose Dr	Resp Rate	Measured respiration rate – capnography (Patient monitor)		
ET Des ET Des End tidal Desflurane anesthetic agent (Patient monitor) ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET N20 End tidal Isoflurane anesthetic agent (Patient monitor) ET N20 End tidal Isoflurane anesthetic agent (Patient monitor) BIS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) FEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Q 2 FLOW Q 3 flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery information of medications Time Drug Name Dose Drug dose Drug dose Drug dose Drug dose Drug dose Unit Drug dose unit Route Drug dose unit Pluid totals Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Volume Fluid name Volume Fluid name Volume Lab results Intraoperative laboratory results Time Time of fluid input or output Fluid Name Fluid name Volume Lab results Lab result description Sample type — arterial/venous Lab result Value Unit Unit Onit of measurement Notes	FiO ₂	Inspired O ₂ (Patient monitor)		
ET ISO End tidal Isoflurane anesthetic agent (Patient monitor) ET NZO End tidal Nitrous oxide (Patient monitor) BIS Bispectral Index (Patient Monitor) SQI Signal quality index of BIS (Patient/BIS monitor) SVO ₂ Mixed venous oxygenation (Patient monitor) SVO ₂ Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) TV Tidal volume (Patient monitor) PIP Peak Inspiratory pressure (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₃ flow rate (Ventilator) N ₄ O FLOW Air flow ate (Ventilator) N ₅ O FLOW Air flow ate (Ventilator) N ₅ O FLOW Air flow ate (Ventilator) N ₅ O FLOW Air flow ate (Ventilator) Delivery time (start / end for infusion medications) Time Delivery time (start / end for infusion medications) Drug Name Dose Drug dose Drug dose Drug dose Drug dose Drug dose unit Drug dose unit Route Drug route Fluid totals Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result Lab result Result value Unit Unit Onit of measurement Notes	ET Sevo	End tidal Sevoflurane anesthetic agent (Patient monitor)		
ET N20 End tidal Nitrous oxide (Patient monitor) BIS Bispectral Index (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) CCC CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PIPP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₃ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₂ O Flow air (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Dose Drug dose Drug dose Drug dose Drug dose Drug dose Unit Drug Ose Drug dose Drug Huid totals Time Time of fluid input or output Fluid Name Fluid totals Time Time of fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type — arterial/venous Lab result Unit Unit of measurement Notes Notes	ET Des	End tidal Desflurane anesthetic agent (Patient monitor)		
BIS Signal quality index of BIS (Patient/BIS monitor) SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SvO2 Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Ayo Flow rate (Ventilator) Air Gardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug Name Drug Name Drug dose Dose Unit Drug dose Dose Unit Drug dose Drug dose Drug dose Drug dose Time Time of fluid input or output Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes	ET ISO	End tidal Isoflurane anesthetic agent (Patient monitor)		
Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP P Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery information of medications Drug Name Drug Mame Dose Drug dose Drug dose Drug dose Drug dose Drug dose Unit Drug dose unit Prud totals Fluid totals Fluid totals Fluid name Fluid name Fluid name Volume Fluid name Fluid Name Fluid name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result Lab result Result value Unit Unit of measurement Notes Notes	ET N20	End tidal Nitrous oxide (Patient monitor)		
SQI Signal quality index of BIS (Patient/BIS monitor) TEMP Temperature (Patient monitor) SVO2 Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP P Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery information of medications Drug Name Drug Name Dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result description Sample type — arterial/venous Lab result Unit Onto measurement Notes Attestations, AIMS Clinical notes, Note option selections	BIS	Bispectral Index (Patient/BIS monitor)		
TEMP Temperature (Patient monitor) SvO ₂ Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₃ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₃ O Flow rate (Ventilator) Time Delivery information of medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Dose Drug Name Dose Drug dose Drug dose Drug dose Drug dose Drug dose Unit Drug route Fluid totals Fluid totals Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result Lab result Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	SQI			
SVO2 Mixed venous oxygenation (Patient monitor) CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₂ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Time Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Aose Dose Unit Drug dose Dose Unit Drug dose Drug dose Drug dose Drug dose Time Time of fluid input or output Fluid totals Fluid totals Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of fluid volume recorded Lab result Lab result description Sample type Sample type a Sample type - arterial/venous Lab result Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	TEMP			
CCO Continuous cardiac output (Patient monitor) TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O.; FLOW O.; flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Time Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections		· · · · · · · · · · · · · · · · · · ·		
TV Tidal volume (Patient monitor) RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₂ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N ₂ O FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₂ O Flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Dose Drug dose Dose Unit Drug dose Dose Unit Drug dose Drug dose Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result description Sample type Sample type - arterial/venous Lab result Result value Unit Unit of measurement Notes	-			
RATE Ventilator rate setting (Ventilator) PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW Air FLOW Air flow rate (Ventilator) N ₂ O FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₂ O Flow rate (Ventilator) Cardiac Rhythm Medications Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug Name Dose Drug dose Drug dose Drug dose Drug route Fluid totals Time Time of fluid input or output Fluid Name Volume Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result Result value Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections		, , ,		
PIP Peak Inspiratory pressure (Ventilator) PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₂ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) Air FLOW N ₂ O FLOW N ₂ O Flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes	• •	,		
PEEP Positive End Expiration Pressure (Ventilator) O ₂ FLOW O ₂ flow rate (Ventilator) Air FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₂ O Flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Time Time of fluid input or output Fluid Name Fluid name Volume Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result description Sample type Sample type – arterial/venous Lab result Unit Ontes Note option selections				
O ₂ FLOW Air FLOW Air flow rate (Ventilator) N ₂ O FLOW N ₂ O FLOW N ₂ O Flow rate (Ventilator) Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug Name Drug Name Drug dose Dose Unit Drug dose Drug dose Drug dose Drug dose Drug route Fluid totals Fluid totals Fluid totals Fluid name Fluid Name Fluid Name Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab result Lab result Sample type Sample type Sample type - arterial/venous Lab result Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections				
Air FLOW N2O FLOW N2O Flow rate (Ventilator) Cardiac Rhythm Medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug dose Drug dose Drug dose unit Route Drug route Fluid totals Time Time of fluid input or output Fluid Name Volume Fluid volume recorded Laboratory results Time Time of taking sample or Lab result time Lab result Sample type Sample type Sample type Sample type Lab result Route Notes Attestations, AIMS Clinical notes, Note option selections		, , ,		
N2O FLOW Cardiac Rhythm Type of cardiac rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Drug Name Drug Name Drug Name Drug Name Drug dose Drug dose Drug dose unit Route Drug route Fluid totals Time Time of fluid input or output Fluid Name Fluid Name Fluid volume recorded Laboratory results Time Time of taking sample or Lab result time Lab result description Sample type Sample type Lab result Lab result Result value Unit Notes Attestations, AIMS Clinical notes, Note option selections		_ ,		
Cardiac Rhythm Medications Delivery information of medications Time Delivery time (start / end for infusion medications) Drug Name Dose Drug Name Dose Dose Unit Route Fluid totals Time Time of fluid input or output Fluid Name Volume Laboratory results Time Time of taking sample or Lab result time Lab result Lab result Sample type Sample type Sample type Lab result Result value Unit Notes		, ,		
MedicationsDelivery information of medicationsTimeDelivery time (start / end for infusion medications)Drug NameDrug NameDoseDrug doseDose UnitDrug dose unitRouteDrug routeFluid totalsFluid totalsTimeTime of fluid input or outputFluid NameFluid nameVolumeFluid volume recordedLaboratory resultsIntraoperative laboratory resultsTimeTime of taking sample or Lab result timeLab resultsLab result descriptionSample typeSample type - arterial/venousLab resultResult valueUnitUnit of measurementNotesAttestations, AIMS Clinical notes, Note option selections				
Time Delivery time (start / end for infusion medications) Drug Name Drug Name Dose Drug dose Dose Unit Drug dose unit Route Drug route Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Drug Name (start / end for infusion medications) Drug Name Drug		· · · · · · · · · · · · · · · · · · ·		
Drug Name Dose Dose Unit Dose Unit Dose Unit Dorug dose Unit Route Drug route Fluid totals Fluid totals Time of fluid input or output Fluid Name Fluid Name Fluid volume recorded Laboratory results Time Time of taking sample or Lab result time Lab results Lab results Sample type Sample type — arterial/venous Lab result Unit Unit of measurement Notes Drug Name Drug Name Drug dose Drug route Drug route Fluid totals Fluid totals Fluid totals Fluid totals Fluid totals Fluid input or output Fluid volume recorded Lab results Fluid sample or Lab result time Lab result description Sample type — arterial/venous Lab result Unit Onit of measurement Notes	Medications	Delivery information of medications		
Dose Unit Drug dose unit Route Drug route Fluid totals Fluid totals Time Time of fluid input or output Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Unit Unit of measurement Notes Notes	Time	Delivery time (start / end for infusion medications)		
Dose Unit Route Drug dose unit Route Fluid totals Fluid totals Fluid name Fluid Name Fluid Name Fluid volume recorded Laboratory results Time Time of taking sample or Lab result time Lab results Lab results Lab results Sample type Sample type Lab result Result value Unit Notes Drug dose unit Drug dose unit Drug dose unit Drug route Fluid volume Fluid totals Fluid input or output Fluid input or output Fluid name Fluid volume recorded Lab results Lab results Lab result time Lab result description Sample type — arterial/venous Lab result Result value Unit of measurement Notes	Drug Name	Drug Name		
Route Drug route Fluid totals Fluid totals Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Time of taking sample or Lab result time Lab results Lab results Lab result description Sample type Sample type — arterial/venous Lab result Unit Unit Of measurement Notes	Dose	Drug dose		
Fluid totals Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Fluid totals Fluid totals Fluid totals Fluid input or output Fluid name Fluid name Sample spections	Dose Unit	Drug dose unit		
Time of fluid input or output Fluid Name Fluid Name Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Time of fluid input or output Fluid name Fluid volume recorded Laboratory results Intraoperative laboratory results Lab result time Lab result value Unit of measurement Notes	Route	Drug route		
Time Time of fluid input or output Fluid Name Fluid name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type – arterial/venous Lab result Result value Unit Unit of measurement Notes Time of taking sample or Lab result time Lab result time Lab result description Sample type Attestations, AIMS Clinical notes, Note option selections	Fluid totals	Fluid totals		
Fluid Name Volume Fluid volume recorded Laboratory results Intraoperative laboratory results Time Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	- 1 - 1 - 1 - 1 - 1 - 1 - 1 - 1 - 1 - 1	Time of fluid input or output		
VolumeFluid volume recordedLaboratory resultsIntraoperative laboratory resultsTimeTime of taking sample or Lab result timeLab resultsLab result descriptionSample typeSample type – arterial/venousLab resultResult valueUnitUnit of measurementNotesAttestations, AIMS Clinical notes, Note option selections	11110	•		
Laboratory results Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type – arterial/venous Lab result Result value Unit Unit of measurement Notes Intraoperative laboratory results Time of taking sample or Lab result time Lab result testing Sample type – arterial/venous Result value Unit of measurement Notes				
Time of taking sample or Lab result time Lab results Lab result description Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections				
Lab results Sample type Sample type — arterial/venous Lab result Result value Unit Unit of measurement Notes Lab result description Attestations, AIMS Clinical notes, Note option selections	Laboratory results	intraoperative laboratory results		
Sample type Sample type – arterial/venous Lab result Result value Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	Time	Time of taking sample or Lab result time		
Lab result Unit Result value Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	Lab results	Lab result description		
Unit Unit of measurement Notes Attestations, AIMS Clinical notes, Note option selections	Sample type	Sample type – arterial/venous		
Notes Attestations, AIMS Clinical notes, Note option selections	Lab result	Result value		
Hotes	Unit	Unit of measurement		
Time Time associated with the note	Notes	Attestations, AIMS Clinical notes, Note option selections		
	Time	Time associated with the note		
Context code Context code associated with a note or its selections	Context code	Context code associated with a note or its selections		
Note content Note description				

Supplementary Table 1. Raw data sources used for machine learning. A wide variety of data sources were used that included text values, time series, discrete values, and quantitative values.

*SpO*₂ desaturation labels

We considered $SpO_2 \le 92\%$ as hypoxemia, which falls between the World Health Organization's recommended intervention level (< 94%) and emergency level (< 90%) (25). Predictions of hypoxemia were made for a window 5 minutes into the future. If the SpO_2 was $\le 92\%$ at any point during those 5 minutes then it was considered a positive label, otherwise it was negative. The machine learning algorithm was trained using these *training labels* on all time points where SpO_2 was not already $\le 92\%$ at that time point.

When evaluating the machine learning algorithm's performance by comparing with anesthesiologists we chose to use a much more stringent definition of hypoxemia. This more stringent definition excludes some of the time points, and leads to a smaller set of *testing labels*. *Testing labels* were positive only if SpO_2 was $\geq 95\%$ for the past 10 minutes and then fell below 92% in the next five minutes (Supplementary Figure 1; left). *Testing labels* were negative only if SpO_2 remained $\geq 95\%$ for the past ten minutes and the next ten minutes (Supplementary Figure 1; right). All the other cases do not have testing labels. This more restrictive labeling scheme ensures that positive testing labels are clear drops in SpO_2 levels, while negative testing labels are clearly not drops in SpO_2 (Supplementary Figure 1).

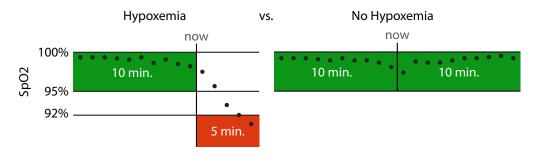

Supplementary Fig 1. Criteria for defining testing labels. When comparing the performance with anesthesiologists, only time points that clearly either desaturate or not are used. Hypoxemia involves dropping from $\geq 95\%$ to $\leq 92\%$ in the next 5 minutes. Not desaturating means remaining $\geq 95\%$ for both the past 10 minutes and the next 10 minutes.

An important point to consider when building labels for health outcome prediction is that anesthesiologist interventions can affect outcomes. It was recently noted by Dyagilev et al. that models can learn when an anesthesiologist is likely to intervene and hence lower the risk of an otherwise high-risk patient (26). This means that patients with low risk (from the model) may still need treatment. To address this, they proposed removing examples from the training set where anesthesiologists have intervened. This allows one to learn a model which predicts patient outcome without intervention. In our case, it is not possible to fully identify when or how an anesthesiologist is intervening (and if that intervention prevented hypoxemia), so we sought to address this issue in two ways:

- 1. It must be recognized that the model predicts hypoxemia when following standard procedures, *not* the occurrence of hypoxemia if the anesthesiologist takes no action to influence hypoxemia. This is a natural assumption in the operating room where interventions that may affect SpO₂ levels are performed frequently.
- 2. By focusing on clear explanations of why a certain risk was predicted we enable anesthesiologists to identify when the algorithm may be basing its risk on their actions vs. when the risk is based on other factors.


Extracting time series features

To make a prediction at an arbitrary point in time, a consistent set of features should be computed that capture the information present in all previous time points. All the data provided about a procedure is associated with a specific date and time. Text data has the time it was provided, minute-by-minute data from the patient monitor has the time at which each measurement was taken, and single point measurements have the times they were recorded.

We summarized these unevenly sampled time registered data into a fixed length feature vector at any point in time using several complementary methods:

- Patient data, procedure information, and pre-operative notes are represented by a "last value" feature, which is zero before any data is recorded and the data's value afterwards.
- Time series data is captured using exponentially decaying weighted average and variance estimates using multiple decay rates. These decay rates specify how much impact each past time point has on the computed mean or variance for the time series. We used 6 second, 1 minute, and 5-minute half-life times to capture both high and low frequency components of the signal in each time series (Supplementary Figure 2).
- Drug dose data is captured using both an exponentially decaying sum, and a time since the last measurement. Decay rates with half-lives of 5 minutes and 60 minutes were used to capture both near term and longer average drug dosing effects.

To ensure that there was enough training data for each feature we removed features that had less than 100 recorded data values for the real-time model, and less than 50 for the initial model. For a full list of the 3,797 features used by Prescience for initial predictions see Supplementary Table 4 (Supp_InitialFeatureTable_4.csv). For the 3,905 features used in intraoperative predictions see Supplementary Table 5 (Supp_RealtimeFeatureTable_5.csv).

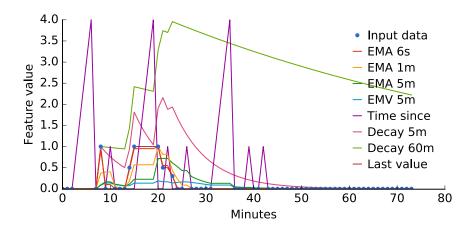

Supplementary Fig 2. Responses of the eight different time series features used in Prescience to a sample set of unevenly reported data values. The blue dots represent the original unevenly sampled data, while curves represent the value of a feature over time. 'EMA' stands for exponential moving average and 'EMV' for exponential moving variance. Both EMA and EMV features are computed over weighted samples, where the weights decay with a specific half-life (6 seconds, 1 minute, or 5 minutes).

Gradient boosting machines for prediction


The features we compute from real-time operating room data have a variety of complex nonlinear interactions. Capturing these requires a model with significant flexibility, and we chose a non-parametric approach called *gradient boosting machines* (19).

We compared the performance of gradient boosting against three baseline methods: Lasso penalized linear logistic regression; a linear SVM autoregressive model previously proposed for predicting hypoxemia based only on the SpO₂ data stream (9); and an unsupervised Parzen window method used previously to predict patient deterioration (17). Gradient boosting machines significantly outperformed all baseline methods for our primary endpoint, real-time hypoxemia prediction (Supplementary Figure 3). For our secondary task of initial prediction gradient boosting machines were only slightly superior (Supplementary Figure 4). The large performance gain of gradient boosting for intraoperative prediction (Supplementary Figure 3) is likely because there are 8 million training samples, while for preoperative predictions (Supplementary Figure 4) there are only 42,000 samples and no time series data. Note that for initial prediction the autoregressive SVM and Parzen window methods were not applicable and hence not evaluated.


Gradient boosting machines are non-parametric models that draw a parallel between boosting and gradient descent in function space. They additively build up simpler models, like boosting, and these models are fit to the gradient of the loss at every data point. The most common type of basic model used is a regression tree because it is both robust to outliers and flexible. Taking some small fraction, η , of many trees fit to the gradient results in many small gradient descent steps in function space.

Fitting the trees is computationally challenging on large datasets so we used XGBoost, a recent high performance implementation of gradient boosting machines (20). For the real-time model we used $\eta = 0.2$ and 1,242 trees, while for the initial model we chose $\eta = 0.1$ and 4,000 trees. Using a smaller η value means more trees are required for fitting, which requires more time to run, but results in a smoother (and generally better) model. For both initial and real-time models we used bagging, where trees were trained on a random 50% subsample of the training data. For the preoperative model the max tree depth was 4 and the minimum child weight of any branch in the trees was 1. For the real-time model the max tree depth was 6 and the minimum child weight of any branch in the trees was 10.

All method parameters were tuned (and methods were chosen) using a validation set of operating room procedures separate from the final test set used for all final performance results. To ensure that there was no bias towards the final test set, the test data was initially compressed and left compressed until after method development was completed.

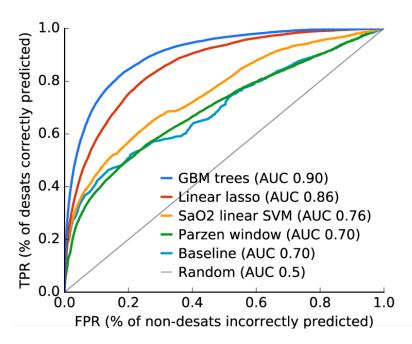

Supplementary Fig 3. Real-time performance of gradient boosting machines vs a linear lasso model, a SVM based on ElMoaqet et al. (9), and an unsupervised Parzen window method used by Tarassenko et al. (17). There are ~8 million training samples and the increased flexibility of gradient boosting trees clearly outperforms the more restrictive linear models.

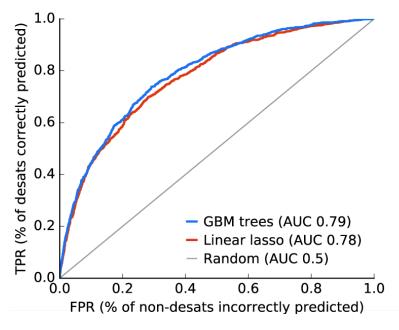

Supplementary Fig 4. Preoperative performance of gradient boosting machines vs a linear lasso model. Given a much smaller preoperative dataset

with ~42,000 training examples the difference between the complex gradient boosting trees model and a linear model becomes small.

Computing feature importance estimates

Understanding why a statistical model has made a specific prediction is a key challenge in machine learning. It engenders appropriate trust in predictions and provides insight into how a model may be improved. However, many complex models with excellent accuracy, such as gradient boosting, make predictions even experts struggle to interpret. This forces a tradeoff between accuracy and interpretability. In response to this we chose to use a model agnostic representation of feature importance, where the impact of each feature on the model is represented using *Shapley values (27)*, theoretically proven unique feature attribution values that have four important properties below. We developed an efficient method to compute the Shapely values (i.e., estimated importance of features to a particular prediction) in a real-time manner.

Shapley values are from the game theory literature and provide a theoretically justified method for allocation of credit among a group (see Eq (1)). In Prescience the group is a set of interpretable model input feature values, and the credit is the value of the prediction made by the model when given those input feature values. Feature impact is defined as the change in prediction probability when a feature is observed vs. unknown. Some feature values have a large impact on the prediction, while others have little impact. The Shapley values $\phi_i(f, x)$, explaining a prediction f(x), are an allocation of credit among the various features (such as age, weight, time series features, etc.), and are the only such allocation that obeys a set of desirable properties (14, 16). Given a prediction function f(x) we can define $f_x(S) = f(x_S)$, where x_S is an input vector equal to x for features in the set S but equal to the missing data value otherwise.

Using this notation, we use the following set of properties:

Efficiency.

$$f(x) = \sum_{i=0}^{M} \phi_i(f, x),$$

where *M* is the number of 'interpretable' inputs, which each correspond to a group of original input features. For Prescience these groups are the sets of features associated with each time series. For instance, the 6 second, 1 minute, and 5 minute moving average features, and the 5-minute moving variance feature from the SpO₂ time series are all considered as a single group. The efficiency assumption forces the model to correctly capture the original predicted value.

Symmetry. If for all subsets S that don't contain i or j

$$f_{x}(S \cup \{i\}) = f_{x}(S \cup \{j\}),$$

then $\phi_i(f, x) = \phi_j(f, x)$. This states that if two features have an identical effect when observed in any situation then the Shapley values for the features must be the same.

Null effects. If for all subsets S that don't contain i

$$f_{\mathcal{X}}(S \cup \{i\}) = f_{\mathcal{X}}(S),$$

then $\phi_i(f, x) = 0$. This states that if a feature has no effect when observed in any situation then its Shapley value must be 0.

Linearity. For any two model prediction functions f and f'

$$\phi_i(f + f', x) = \phi_i(f, x) + \phi_i(f', x).$$

This states that the effect a feature has on the sum of two functions is the effect it has on one function plus the effect it has on the other.

Shapley showed that *only one allocation of credit satisfies these properties* and that allocation is the one given by the Shapley values.

Given a specific prediction f(x) we can compute the Shapley values using a weighted sum:

$$\phi_{i}(f,x) = \sum_{S \subseteq S_{all} \setminus \{i\}} \frac{|S|! (M - |S| - 1)!}{M!} [f_{x}(S \cup \{i\}) - f_{x}(S)]$$

$$= \sum_{S \subseteq S_{all} \setminus \{i\}} \frac{1}{(M \ choose \ |S|)(M - |S|)} [f_{x}(S \cup \{i\}) - f_{x}(S)].$$
Eq (1)

In practice, there are far too many terms to evaluate this sum completely, so we instead calculate it by only selecting the highest weight terms. This selection can be done greedily, stopping once the variance of the estimate becomes small enough.

To compute the Shapley values of each prediction we need to estimate the predictions of the model when specific input features are missing. Since the model was not trained to support missing values we approximate what the model would predict (if retrained on that subset of input features) by sampling from the training data set and replacing the missing features with the values they would have had in that sample. By averaging many such samples, we can estimate the expected value of $f_x(S)$ only using evaluations of $f_x(S_{all})$ where no features are missing.

The approach above requires nested sampling, once to estimate the Shapley value and then from each sample we again sample to estimate $f_x(S)$ and $f_x(S \cup \{i\})$. To reduce the number of samples in the inner step, we used k-medians to generate 20 medians of the entire dataset, and then performed a weighted evaluation for only these 20 summary inputs as an approximation for the entire dataset. This removes the need for nested sampling.


One final extension of Shapley value estimation we found to be helpful in Prescience was a non-linear link function *h* such that:


$$h(f(x)) = \sum_{i=0}^{M} \phi_i(f, x).$$

Since Prescience uses logistic regression the use of a h = logit link function transforms the output space from probabilities to log odds. Assuming the importance of features is additive in the log odds space is much more natural than assuming they are additive in the space of probabilities (which must fall between 0 and 1). The same reasoning also drives the use of the logit link function during standard logistic regression.

We were able to get stable feature importance estimates for thousands of features in less than 5 seconds on our server. We compared these theoretically grounded explanations with a simple estimate of feature

importance to verify they showed reasonable consistency. The simple method we chose was to replace a single feature group with random values from other samples in the data set, and determine the average model output over different possible samplings. We then subtracted this mean value from the original model prediction to get a difference from a prediction with a typical value of that feature vs. the current value. This simple method is not very scalable and does not account for interactions with other features, but is useful to compare with the Prescience explanations to ensure the Prescience estimates of feature effects are consistent with an intuition of how much a feature's change from its typical value effects the current risk of hypoxemia (Supplementary Figure 5).

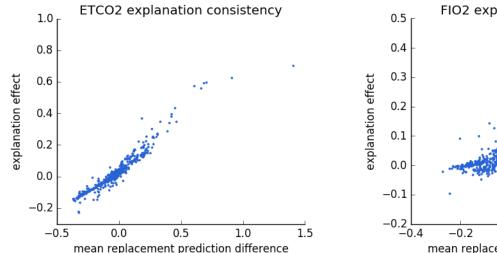

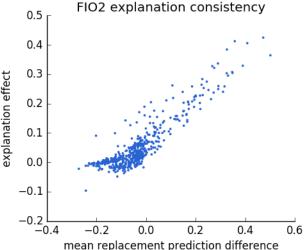

Supplementary Fig 5. Consistency of Prescience explanation effects. Comparing the Prescience explanation effects with the difference between the current model output and the outputs when a specific feature is replaced with its typical value. The effects are from the cases shown to anesthesiologists during testing in Figure 2. The strong correlation ($R^2 = 0.92$ for ETCO₂ and $R^2 = 0.81$ for FIO₂) demonstrates the consistency of Prescience explanation effect sizes with the intuitive notion of the change in model out from replacing a feature's value with a typical value. Note that both values are shown explaining the additive portion of the classification model (inside the logistic function).

Physician evaluation


The potential benefit Prescience provides to physicians was evaluated using previously recorded procedures. Both before a procedure begins, and at several time points during the operation all the available electronically recorded data was shown to the anesthesiologist and they were asked to predict if a desaturation (as defined above) will occur in the next 5 minutes (Supplementary Figures 6, 7, and 8). For half of the procedures anesthesiologists are given Prescience explained risks, and for the other half they are given the same data, but without any Prescience assistance. In both cases anesthesiologists are asked to provide a fold change in the risk that desaturation will occur.

The test procedures were divided into two equal sized groups, replicate 1 and replicate 2. Anesthesiologists were also divided into two groups, A and B. Group A was given Prescience assistance on replicate 1 but not on replicate 2, while group B was given Prescience assistance on replicate 2 but not replicate 1. After randomly assigning anesthesiologists to groups, three anesthesiologists from group A completed the evaluation and two anesthesiologists from group B. We pooled the results within each


group and between groups, and the results of this evaluation are shown in Figure 2. The order in which anesthesiologists were presented with cases was random across both replicate sets.



Supplementary Fig 6. Sample of a physician's test interface for preoperative prediction of risk. Prescience assistance is given for this preoperative prediction. Anesthesiologists choose an estimated relative risk by moving the given slider, then recording the score.

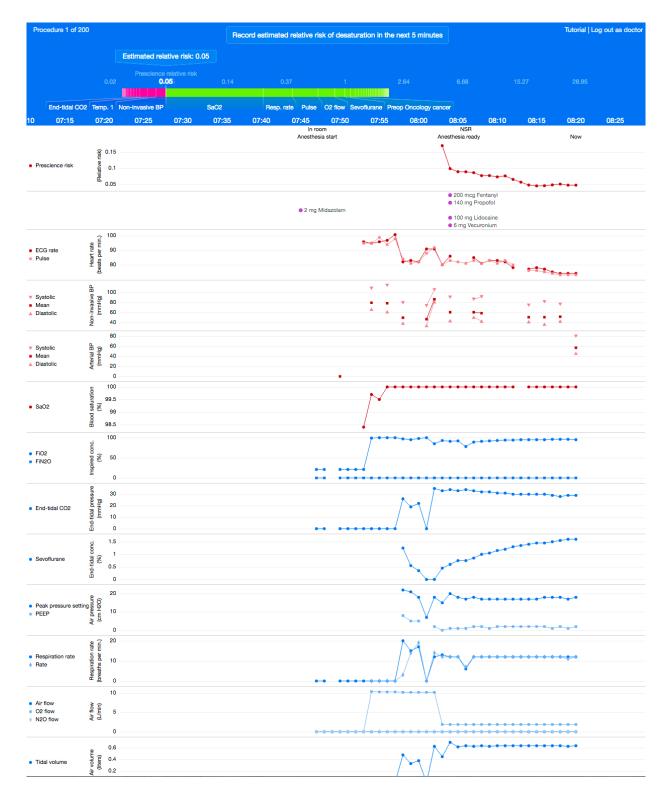

Supplementary Fig 7. Sample of a physician's test interface for real-time prediction of risk. Prescience assistance is given for this intraoperative prediction. Anesthesiologists choose an estimated relative risk by moving the given slider. The preoperative data is also shown lower down in the interface just as illustrated in Supplementary Figure 6.

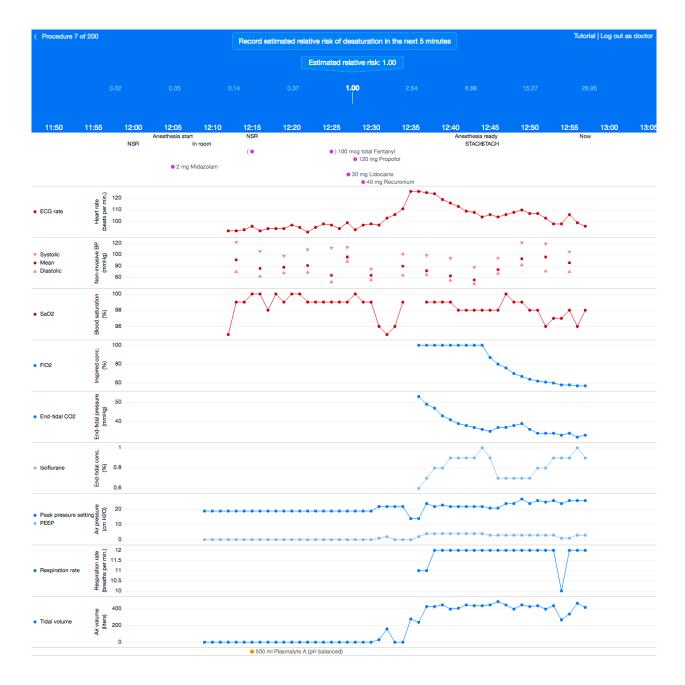

Supplementary Fig 8. Sample of a physician's test interface for real-time prediction of risk without Prescience assistance. This is similar to Supplementary Figure 7 but this time the anesthesiologist must choose a risk on their own using only the original data without help from Prescience.

	Initial predictions					
	Anesthesiolog ist 1	Anesthesiolo gist 2	Anesthesiologi st 3	Anesthesiologist 4	Doctor composite	
1	surgical procedure	BMI	BMI (and by default height and weight)	Baseline saturation (patient on home O ₂ /room air)	BMI	
2	BMI	Weight	History of OSA	Pulmonary disease (Fibrosis, asthma, COPD - Chronic Obstructive Pulmonary Disease)	Lung disease	
3	age	Smoker	Type of surgery	Obesity or high BMI (Body mass index)	ASA code	
4	pre-existing heart or lung disease	Presence of lung disease	History of Asthma	Certain planned procedures (Thoracic surgery)	Asthma	
5	ASA status	History of difficult airway	Functional Status	Anesthesia method (sedation)	Age	
6				Children	Anesthesia type	
7				Pregnancy	CPT code	
8				High metabolic status (sepsis)		

Supplementary Table 2. Physician's responses listing the most important factors when assessing the risk for hypoxemia at any point during an upcoming procedure. Four anesthesiologists were asked to list the most important factors to consider for hypoxemia. The responses of all 4 anesthesiologists were then compiled into a single composite list. For a comparison with features chosen by Prescience, see Figure 4 in the main text.

	Real-time predictions						
	Anesthesiolo gist 1	Anesthesiologist 2	Anesthesiologis t 3	Anesthesiologist 4	Anesthesiologis t composite		
1	surgical procedure	Trend of SPO ₂ if downward trend more likely	Tidal Volume	Trends in SpO ₂	SpO ₂		
2	BMI	Tidal volume - if trending down suggests a problem.	Respiratory Rate	Surgical interventions (1-lung ventilation, insufflation of the abdomen, pushing the lung)	Tidal Volume (TV)		
3	FiO ₂	Increasing Peak inspiratory pressure	SpO_2	High airway pressure	Peak Inspiratory Pressure (PIP)		
4	minute ventilation/tid al volume	Increasing Co ₂	End Tidal CO ₂	Surgical and Anesthesia events (bronchospasm, pneumothorax, pulmonary embolism, mucus plug, airway trapping, fighting against vent etc.)	End tidal CO ₂		
5	recent anesthetic drug administration	Increasing heart rate	Peak Inspiratory Pressure (PIP)	Patient position (head down)	BMI		
6			Tachycardia	Anesthesia machine malfunction (tube kinking)	Resp rate		
7			BMI (and by default height and weight)		FiO2		
8			History of OSA		Heart rate		
9			Type of surgery				
10			History of Asthma				

Supplementary Table 3. Physician's responses listing the most important factors when assessing the risk for hypoxemia in the next 5 minutes. Four anesthesiologists were asked to list the most important factors to consider for hypoxemia. The responses of all 4 anesthesiologists were then compiled into a single composite list. For a comparison with features chosen by Prescience, see Figure 4 in the main text.

Supplementary Materials

Supplementary Fig 1. Criteria for defining testing labels. When comparing the performance with anesthesiologists, only time points that clearly either desaturate or not are used. Hypoxemia involves dropping from $\geq 95\%$ to $\leq 92\%$ in the next 5 minutes. Not desaturating means remaining $\geq 95\%$ for both the past 10 minutes and the next 10 minutes.

Supplementary Fig 2. Responses of the eight different time series features used in Prescience to a sample set of unevenly reported data values. The blue dots represent the original unevenly sampled data, while curves represent the value of a feature over time. 'EMA' stands for exponential moving average and 'EMV' for exponential moving variance. Both EMA and EMV features are computed over weighted samples, where the weights decay with a specific half-life (6 seconds, 1 minute, or 5 minutes).

Supplementary Fig 3. Real-time performance of gradient boosting machines vs a linear lasso model, a SVM based on ElMoaqet et al. (9), and an unsupervised Parzen window method used by Tarassenko et al. (17). There are \sim 8 million training samples and the increased flexibility of gradient boosting trees clearly outperforms the more restrictive linear models.

Supplementary Fig 4. Initial performance of gradient boosting machines vs a linear lasso model. Given a much smaller preoperative dataset with ~42,000 training examples the difference between the complex gradient boosting trees model and a linear model becomes small.

Supplementary Fig 5. Consistency of Prescience explanation effects. Comparing the Prescience explanation effects with the difference between the current model output and the outputs when a specific feature is replaced with its typical value. The effects are from the cases shown to anesthesiologists during testing in Figure 2. The strong correlation ($R^2 = 0.92$ for ETCO₂ and $R^2 = 0.81$ for FIO₂) demonstrates the consistency of Prescience explanation effect sizes with the intuitive notion of the change in model out from replacing a feature's value with a typical value. Note that both values are shown explaining the additive portion of the classification model (inside the logistic function).

Supplementary Fig 6. Sample of a physician's test interface for initial prediction of risk. Prescience assistance is given for this preoperative prediction. Anesthesiologists choose an estimated relative risk by moving the given slider, then recording the score.

Supplementary Fig 7. Sample of a physician's test interface for real-time prediction of risk. Prescience assistance is given for this intraoperative prediction. Anesthesiologists choose an estimated relative risk by moving the given slider. The preoperative data is also shown lower down in the interface just as illustrated in Supplementary Figure 5.

Supplementary Fig 8. Sample of a physician's test interface for real-time prediction of risk without Prescience assistance. This is similar to Supplementary Figure 6 but this time the anesthesiologist must choose a risk on their own using only the original data without help from Prescience.

Supplementary Table 1. Raw data sources used for machine learning. A wide variety of data sources were used that included text values, time series, discrete values, and quantitative values.

Supplementary Table 2. Physician's responses listing the most important factors when assessing the risk for hypoxemia at any point during an upcoming procedure. Four anesthesiologists were asked to list the most important factors to consider for hypoxemia. The responses of all 4 anesthesiologists were then compiled into a single composite list. For a comparison with features chosen by Prescience, see Figure 4 in the main text.

Supplementary Table 3. Physician's responses listing the most important factors when assessing the risk for hypoxemia in the next 5 minutes. Four anesthesiologists were asked to list the most important factors to consider for hypoxemia. The responses of all 4 anesthesiologists were then

compiled into a single composite list. For a comparison with features chosen by Prescience, see Figure 4 in the main text.

Supplementary Table 4. Preoperative features used by Prescience. An enumeration of all the 3,797 features used for preoperative predictions (Supp_InitialFeatureTable_4.csv).

Supplementary Table 5. Intraoperative features used by Prescience. An enumeration of all the 3,905 features used for intraoperative predictions (Supp_RealtimeFeatureTable_5.csv).

Code availability

All modeling and processing code is available from the authors upon request. However, note that patient privacy prevents training data from accompanying the code.

Data availability

The operating room datasets from participating hospitals are not publically available due to patient privacy concerns. For more information on obtaining additional summaries of the data please contact the corresponding author.

Funding

This work was supported by a National Science Foundation (NSF) DBI-135589, NSF Graduate Research Fellowship, and a UW eScience/ITHS seed grant *Machine Learning in Operating Rooms* (06-1019).

References and Notes:

- 1. A. B. Haynes, T. G. Weiser, W. R. Berry, S. R. Lipsitz, A.-H. S. Breizat, E. P. Dellinger, T. Herbosa, S. Joseph, P. L. Kibatala, M. C. M. Lapitan, A. F. Merry, K. Moorthy, R. K. Reznick, B. Taylor, A. A. Gawande, Safe Surgery Saves Lives Study Group, A Surgical Safety Checklist to Reduce Morbidity and Mortality in a Global Population, *N. Engl. J. Med.* **360**, 491–499 (2009).
- 2. A. A. Gawande, E. J. Thomas, M. J. Zinner, T. A. Brennan, The incidence and nature of surgical adverse events in Colorado and Utah in 1992, *Surgery* **126**, 66–75 (1999).
- 3. A. K. Kable, R. W. Gibberd, A. D. Spigelman, Adverse events in surgical patients in Australia., *Int. J. Qual. Heal. care J. Int. Soc. Qual. Heal. Care* **14**, 269–76 (2002).
- 4. B. G. Nair, E. Gabel, I. Hofer, H. A. Schwid, M. Cannesson, Intraoperative Clinical Decision Support for Anesthesia, *Anesth. Analg.* **124**, 603–617 (2017).
- 5. G. D. H. & P. J. Lena Maier-Hein, Swaroop S. Vedula, Stefanie Speidel, Nassir Navab, Ron Kikinis, Adrian Park, Matthias Eisenmann, Hubertus Feussner, Germain Forestier, Stamatia Giannarou, Makoto Hashizume, Darko Katic, Hannes Kenngott, Michael Kranzfelder, Anand Malpani, Surgical data science for next-generation interventions, *Nat. Biomed. Eng.* **1**, 691–696 (2017).
- 6. C. M. Dunham, B. M. Hileman, A. E. Hutchinson, E. A. Chance, G. S. Huang, Perioperative hypoxemia is common with horizontal positioning during general anesthesia and is associated with major adverse outcomes: a retrospective study of consecutive patients, *BMC Anesthesiol.* **14**, 43 (2014).
- 7. L. Strachan, D. W. Noble, Hypoxia and surgical patients--prevention and treatment of an unnecessary cause of morbidity and mortality., *J. R. Coll. Surg. Edinb.* **46**, 297–302 (2001).
- 8. J. M. Ehrenfeld, L. M. Funk, J. Van Schalkwyk, A. F. Merry, W. S. Sandberg, A. Gawande, The incidence of hypoxemia during surgery: evidence from two institutions, *Can. J. Anesth. Can. d'anesthésie* **57**, 888–897 (2010).

- 9. H. ElMoaqet, D. M. Tilbury, S. K. Ramachandran, Multi-Step Ahead Predictions for Critical Levels in Physiological Time Series, *IEEE Trans. Cybern.* **46**, 1704–1714 (2016).
- 10. Z. C. Lipton, D. C. Kale, R. C. Wetzell, Phenotyping of Clinical Time Series with LSTM Recurrent Neural Networks, (2015) (available at http://arxiv.org/abs/1510.07641).
- 11. K. E. Henry, D. N. Hager, P. J. Pronovost, S. Saria, A targeted real-time early warning score (TREWScore) for septic shock, *Sci. Transl. Med.* **7** (2015).
- 12. S. Saria, A. K. Rajani, J. Gould, D. Koller, A. A. Penn, Integration of Early Physiological Responses Predicts Later Illness Severity in Preterm Infants, *Sci. Transl. Med.* **2**, 48ra65-48ra65 (2010).
- 13. R. C. Deo, Machine Learning in Medicine, *Circulation* 132 (2015).
- 14. E. Štrumbelj, I. Kononenko, Explaining prediction models and individual predictions with feature contributions, *Knowl. Inf. Syst.* **41**, 647–665 (2014).
- 15. M. T. Ribeiro, S. Singh, C. Guestrin, Why Should I Trust You? Explaining the Predictions of Any Classifier, doi:10.1145/2939672.2939778.
- 16. S. Lundberg, S.-I. Lee, An unexpected unity among methods for interpreting model predictions, (2016) (available at https://arxiv.org/abs/1611.07478).
- 17. L. Tarassenko, A. Hann, D. Young, Integrated monitoring and analysis for early warning of patient deterioration, *Br. J. Anaesth.* **97**, 64–68 (2006).
- 18. American Medical Association, Current procedural terminology: CPT (2007).
- 19. J. H. Friedman, Greedy function approximation: A gradient boosting machine., *Ann. Stat.* **29**, 1189–1232 (2001).
- 20. T. Chen, C. Guestrin, XGBoost: A Scalable Tree Boosting System, doi:10.1145/2939672.2939785.
- 21. F. Lumachi, B. Marzano, G. Fanti, S. M. M. Basso, F. Mazza, G. B. Chiara, Relationship between body mass index, age and hypoxemia in patients with extremely severe obesity undergoing bariatric surgery., *In Vivo* **24**, 775–7.
- 22. S. M. Kendale, J. D. Blitz, Increasing body mass index and the incidence of intraoperative hypoxemia, *J. Clin. Anesth.* **33**, 97–104 (2016).
- 23. Defining Adult Overweight and Obesity | Overweight & DC (available at https://www.cdc.gov/obesity/adult/defining.html).
- 24. J. Guay, E. A. Ochroch, in *Cochrane Database of Systematic Reviews*, J. Guay, Ed. (John Wiley & Sons, Ltd, Chichester, UK, 2015), p. CD011151.
- 25. World Health Organization, *Pulse Oximetry Training Manual* (2011).
- 26. K. Dyagilev, S. Saria, Learning (predictive) risk scores in the presence of censoring due to interventions, *Mach. Learn.* **102**, 323–348 (2016).
- 27. A. E. Roth, The Shapley Value Cambridge University Press*Cambridge Univ. Press* (1988) (available at http://www.cambridge.org/catalogue/catalogue.asp?isbn=0511829728).
- 28. W. S. Cleveland, R. Mcgill, Graphical Perception: Theory, Experimentation, and Application to the Development of Graphical Methods, *Source J. Am. Stat. Assoc.* **79**, 531–554 (1984).